

TORS

Orthopaedic Revision System

ARCHIVED COPY.

User Manual

DO NOT USE



CONTENTS

CONVENTIONS USED IN THIS DOCUMENT 65

Warning:..... 65

Caution:..... 65

CONTACT INFORMATION 65

INTENDED USE 76

PLEASE NOTE:..... 76

INDICATIONS..... 76

CONTRAINDICATIONS..... 76


THE TORS SYSTEM 87

GENERATOR..... 87

Front Panel

Rear Panel.....

TORS Transduce

SAFETY .....

WARNINGS.....

CAUTIONS.....

COMPLICATIONS AND PI

INSTRUCTIONS FOR USE..... 1514

SETTING UP THE TORS SYSTEM 1514

Power up the Generator..... 1514

Attach the “Cement” footswitch to the Rear Panel..... 1514

Optional - Attach the “Soft Tissue” footswitch to the Rear Panel..... 1514

Assemble the Cement Transducer + Probe (Applied Part)..... 1615

Optional - Attach the Soft Tissue Transducer (Applied Part)..... 1615

Cement Transducer - Connect 1716

Cement Transducer - Initialise..... 1716

Cement Transducer - Activate..... 1716

Cement Removal – Irrigation..... 1817


Fume Extraction..... 1817

Soft Tissue Transducer - Connect 1918

ARCHIVED COPY:
DO NOT USE

<i>Soft Tissue Transducer – Initialize</i>	1918
<i>Soft Tissue Transducer – Change Power Level</i>	1918
<i>Soft Tissue Transducer – Activate</i>	2019
<i>Turn Off TORS</i>	2019
FAULT / WARNING INDICATORS: (SEE ALSO APPENDIX 4: TONES AND BEEPS)	2120
DECONTAMINATION	2523
REPROCESSING GUIDELINES: IMMEDIATELY AFTER USE	2724
<i>After Use</i>	2724
<i>Single Use Items</i>	2724
MANUAL CLEANING	2724
AUTOMATED CLEANING	272424
CLEANING INSPECTION	2724
WRAPPING	282425
STERILIZATION	2825
<i>End of Life Management</i>	2825
<i>Opened in Error</i>	2825
<i>Cleaning the Generator</i>	3026
<i>Cleaning the Foc</i>	
EQUIPMENT CARE	
<i>Probe Inspection</i>	
<i>Transducer / Pro</i>	
<i>Storage of Equip.</i>	
<i>Maintenance</i>	
<i>Technical Suppor</i>	
APPENDICES	3228
APPENDIX 1: MARKINGS ON THE TORS FRONT AND BACK PANELS	3228
APPENDIX 2: TORS PARTS LIST	35313130
APPENDIX 3: ELECTROMAGNETIC INTERFERENCE	36323231
APPENDIX 4: TONES AND BEEPS	39353534
APPENDIX 5: REAR PANEL MESSAGES	40363635
APPENDIX 6: TECHNICAL SPECIFICATION	41373736

**ARCHIVED COPY:
DO NOT USE**

APPENDIX 7: WARRANTY STATEMENT	42383837
CONVENTIONS USED IN THIS DOCUMENT	4
<i>Warning:</i>	4
<i>Caution:</i>	4
CONTACT INFORMATION	4
INTENDED USE	5
PLEASE NOTE:	5
INDICATIONS	5
CONTRAINDICATIONS	5
THE TORS SYSTEM	6
GENERATOR	6
<i>Front Panel</i>	
<i>Rear Panel</i>	
<i>TORS Transducer</i>	
SAFETY 	
WARNINGS	
Cautions — 10	
INSTRUCTIONS FOR	
SETTING UP THE TORS SYSTEM	11
<i>Power up the Generator</i>	11
<i>Attach the “Cement” footswitch to the Rear Panel</i>	11
<i>Optional – Attach the “Soft Tissue” footswitch to the Rear Panel</i>	11
<i>Assemble the Cement Transducer + Probe</i>	12
<i>Optional – Attach the Soft Tissue + Single Use Handpiece</i>	12
<i>Cement Transducer – Connect</i>	13
<i>Cement Transducer – Initialise</i>	13
<i>Cement Transducer – Activate</i>	13
<i>Cement Removal – Irrigation</i>	14
<i>Soft Tissue Transducer – Connect</i>	15
<i>Soft Tissue Transducer – Initialise</i>	15
<i>Soft Tissue Transducer – Change Power Level</i>	15

**ARCHIVED COPY:
DO NOT USE**

Formatted: Normal, Indent: First line: 0.42 cm

<i>Soft Tissue Transducer—Activate</i>	16
<i>Turn Off TORS</i>	16
FAULT / WARNING INDICATORS: (SEE ALSO APPENDIX 4: TONES AND BEEPS).....	17
DECONTAMINATION	20
REPROCESSING GUIDELINES: IMMEDIATELY AFTER USE.....	21
<i>After Use</i>	21
MANUAL CLEANING.....	21
AUTOMATED CLEANING.....	21
CLEANING INSPECTION.....	21
WRAPPING.....	22
STERILIZATION.....	22
<i>Opened in Error</i>	22
<i>Cleaning the Generator</i>	23
<i>Cleaning the Footswitch</i>	23
EQUIPMENT CARE	24
<i>Probe Inspection</i>	
<i>Transducer / Pro</i>	
<i>Storage of Equip.</i>	
<i>Maintenance</i>	
<i>Technical Support</i>	
APPENDICES	
APPENDIX 1: MARK.....	
APPENDIX 2: TORS PARTS LIST.....	27
APPENDIX 3: ELECTROMAGNETIC INTERFERENCE.....	28
APPENDIX 4: TONES AND BEEPS.....	31
APPENDIX 5: REAR PANEL MESSAGES.....	32
APPENDIX 6: TECHNICAL SPECIFICATION.....	33
APPENDIX 7: WARRANTY STATEMENT.....	34

**ARCHIVED COPY:
DO NOT USE**

Conventions used in this document

Warning:



A statement that if not strictly followed could result in injury or loss of life, or a statement that is otherwise required to be a warning.

Caution:

A statement that if not strictly followed could result in damage to the equipment

Contact Information

If the unit requires repa

RSL Ltd.
Bremridge House
Ashburton
S. Devon
TQ13 7JX
UK

**ARCHIVED COPY:
DO NOT USE**

Intended use

TORS is an ultrasonic surgical device intended to be used for removal of polymethylmethacrylate (PMMA) bone cement in orthopaedic applications. It is also intended to be used for cutting and cauterising soft tissue during skin and muscle incision to access bone joints in limbs.

Please note:

- This document
- The safety & the surgeon

Indications

TORS is indicated for orthopaedic control and r

Contraindications

- Do not use, if in the judgment of the surgeon, the use of ultrasonic surgical techniques are not in the best interest of the patient.
- Do not use for incising bone.
- The Soft Tissue Transducer is not indicated for the cutting of vessels exceeding 1.5mm in diameter or bone tissue.

Please read all Warning & Cautions contained in this document

**ARCHIVED COPY:
DO NOT USE**

The TORS System

Generator

Front Panel

The screen is highly visible, especially in a darkened operating theatre, indicates operating mode selected and gives concise instructions on action required in the event of disruption to normal operation.



Rear Panel

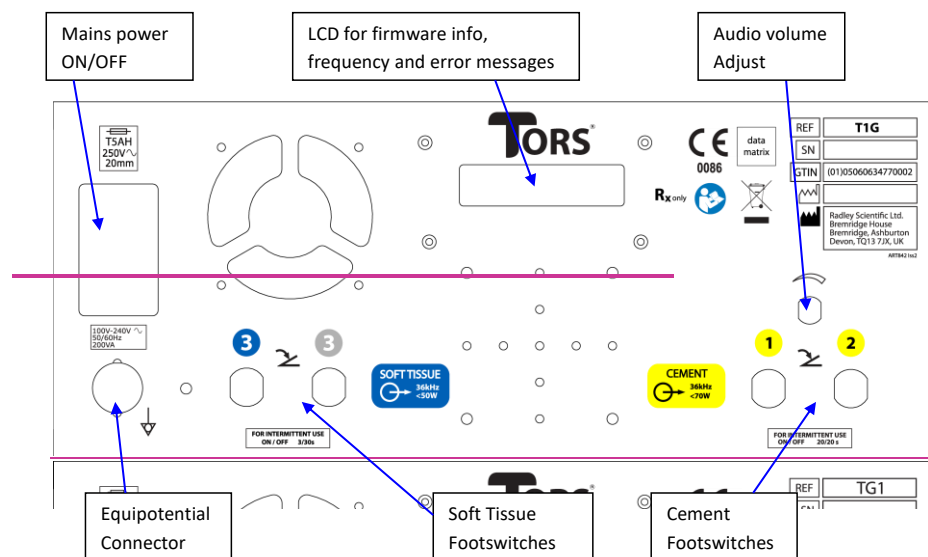


Figure 2: TORS Generator Rear Panel

There is only one user control on the generatorGenerator itself:

- Volume control (on the rear panel)

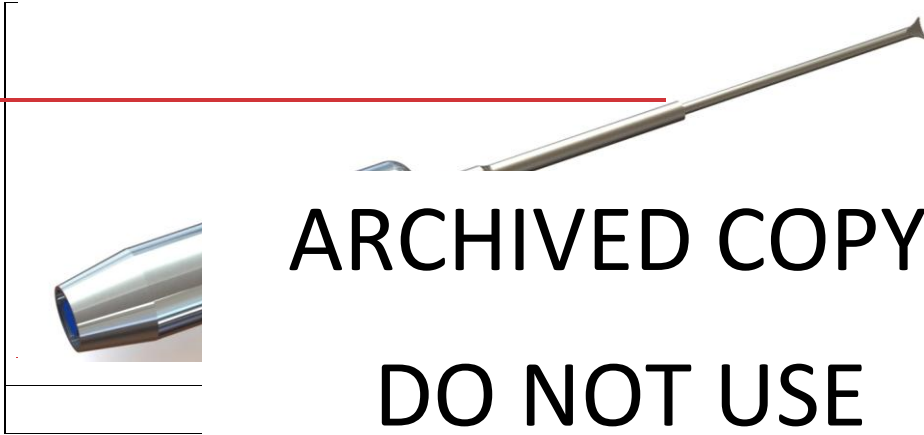


Figure 4: TORS Soft Tissue Transducer + Probe

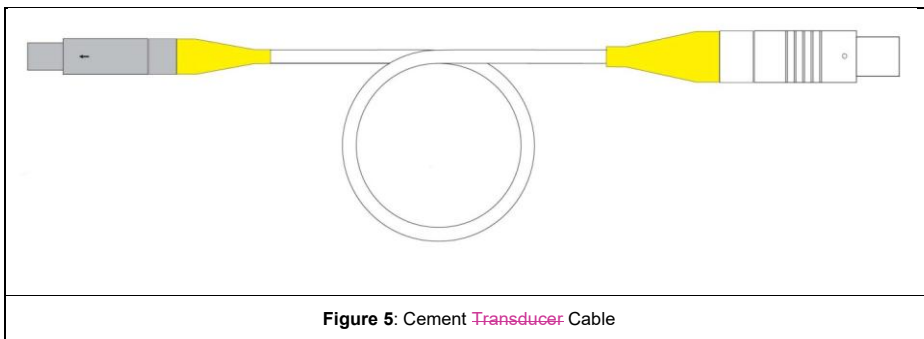


Figure 5: Cement Transducer Cable



**ARCHIVED COPY:
DO NOT USE**



Figure 7 – Axial Grip with TDCR and Probe fitted

Safety

Warnings

- This device should only be used by surgeons who are (1) trained in the types of surgical procedures that are to be carried out and (2) trained in the specific use of ultrasonic surgical instruments.
- Do not operate TORS in a potentially explosive or flammable area, or in oxygen rich environments. Note that if the Probe touches anything metal while activated, sparks may be produced.
- The use of TORS in cases where a patient or an operator of the unit has been fitted with a cardiac pacemaker is left to the discretion of the consultant in charge of the procedure although no electrical flow is made through or into the patient and many clinical papers on the use of ultrasound in the proximity of cardiac pacemakers have been published.
- Appropriate protective measures, including smoke evacuation, should be taken to protect users from any smoke, or other aerosols, produced by the use of surgical ultrasound.
- Use extreme caution when using in the vicinity of nerves.
- Handle the acoustics (Transducers and Probes) with care.
- Do not attempt to modify the acoustics. No modification of this equipment is allowed.
- The use of transducers or handpieces not supplied as part of the TORS system may damage the ~~generator~~Generator and create a safety hazard for the operator and patient.
- Avoid touching or holding the blade at the end of the energized Probe. It is designed to produce heat!
- Other than the activated tip, do not allow probe to be in contact with tissue.
- Irrigate the femoral canal after each activation.
- When guiding the TORS Cement Probe through PMMA, care must be taken not to exert excessive force in order to maintain control of the direction of travel of the Probe.
- Use of excessive force in vicinity of thin bone may cause perforation.
- Avoid resting the blade on skin or other tissue for at least 10s after cessation of energizing as it will become hot while cutting. Ultrasonic systems dissipate heat quickly. (User may quench hot Probe tip / blade in saline)
- Care should be taken when in contact with tissue between activations, in case accidental activation should occur.
- If the Probe comes into contact with bone, using recommended cooling strategies a bone temperature ~~rise~~ rise of up to 5¹°C may be seen.
- Allowing active tip to contact bone may cause damage to the bone.
- Wherever possible avoid contact between the side of active Probe and patient tissue.
- Mains ~~i~~solation is achieved by use of the double pole switch located on the rear panel. DO NOT position the equipment in such a way as to make access to this disconnection switch difficult.
- To avoid the risk of electric shock, this equipment must only be connected to a mains supply with protective earth.
- TORS is not MR safe or MR compatible.
- Do not re-sterilize or reuse any Single Use parts.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of TORS, including cables specified by RSL. Otherwise, degradation of the performance of this equipment could result.
- TORS should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary TORS should be observed to verify normal operation in the configuration in which it will be used
- The use of accessories, transducers and cables other than those specified may result in increased emissions, decreased immunity or improper operation.
- There are no user serviceable parts.

Complications and Potential Side-effects

General risks and complications of arthroplasty surgery

- Bone loss during prosthesis and/or cement removal
- Cortical perforation
- Fracture in the bone around the artificial joint
- Recovery time, pain and potential for arthritis associated with deliberate osteotomy
- Nerve/vascular lesions
- Rebleeding or haematoma / seroma
- Ligament, artery or nerve damage in the area around the knee joint
- Blood clots or deep vein thrombosis or pulmonary embolism
- Haemorrhage caused by tissue injury or ineffective haemostasis
- Bone/joint infection
- Late infection
- Infection of the wound and/or the tissue around the artificial joint
- Thromboembolism
- Implants or metal parts remaining in joint
- Fatigue fracture
- Excess bone forming around the artificial joint (knee) restricting movement
- Excess scar tissue forming and restricting movement
- Kneecap or hip dislocation
- Numbness in the vicinity of the wound scar
- Loosening of the artificial joint (normally after 10-15 years)
- Joint wear and tear
- Joint stiffening
- Bone cement implantation syndrome
- Hypoxia
- Hypotension
- Cardiac arrhythmias
- Cardiovascular collapse
- Thermal injury due to cement polymerization (setting) reaction
- Contact inhalational exposure to PMMA monomer leading to hypersensitivity
- Asthmatic reactions
- Neurological symptoms
- Localised irritation
- Personal risks from anaesthetic or the surgery itself

Potential side-effects of ultrasonic arthroplasty surgery

- Necrosis and clinical dysfunction in femoral cortex and surrounding soft tissue due to heat transmission from ultrasonic probe/cement interface
- Intramedullary canal perforation

Formatted: Font: (Default) Arial, 10 pt

Formatted: Font: (Default) Arial, 10 pt

Formatted: Font: 10 pt

Formatted: Font: 10 pt

Formatted: Indent: Left: 0.63 cm, No bullets or numbering

Formatted: Line spacing: 1.5 lines

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Formatted: Line spacing: 1.5 lines

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Formatted: Line spacing: 1.5 lines

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Formatted: Line spacing: 1.5 lines

Formatted: Indent: Left: 1.9 cm, Line spacing: 1.5 lines, No bullets or numbering

Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering

Formatted: Line spacing: 1.5 lines

Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Formatted: Font: 10 pt

Formatted: Font: 10 pt

- Different tissues and bone cement will be heated to different degrees with the same ultrasound parameters
- Each cell type has a different susceptibility to thermal injury
- Pathologic humeral fracture
- Microscopic bone cracks
- Plume production (95% water, 5% cell debris)
- Radial nerve palsy
- Risk of cross-infection from resterilized reusable parts of the system

- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt

General complications of powered dissection in open soft tissue surgery

- Acute or delayed bleeding (blood loss), organ perforation or fistula formation
- Haemorrhage caused by tissue injury or ineffective haemostasis
- Necrosis due to excessive heat transmission from the dissector
- Infections at the site of application
- Nerve lesions/damage

- Formatted: Font: 10 pt
- Formatted: Indent: Left: 1.9 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Line spacing: 1.5 lines
- Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Line spacing: 1.5 lines
- Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Potential side-effects of ultrasonic powered dissection in open soft tissue surgery

- Necrosis in soft tissue due to cavitation/shock waves or excessive heat from ultrasonic probe
- Visibility reduction due to plume production (95% water, 5% cell debris) or smoke production
- Risk of transmission of mutagenic or toxic chemicals in the plume
- Risk of cross-infection from resterilized reusable parts of the system
- Device breakage and loss of components in the body cavity

- Formatted: Font: (Default) Arial, 9 pt
- Formatted: Normal, Indent: Left: 1.27 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Indent: Left: 0.75 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Normal, Left, Indent: Left: 0.75 cm, Line spacing: 1.5 lines
- Formatted: Font: 10 pt
- Formatted: Line spacing: 1.5 lines
- Formatted: Font: 9 pt
- Formatted: Font: 9 pt
- Formatted: Font: (Default) Arial, 9 pt

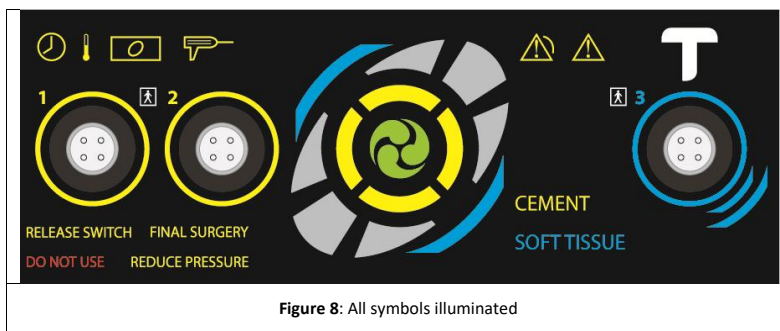
Instructions for Use

Setting Up the TORS System

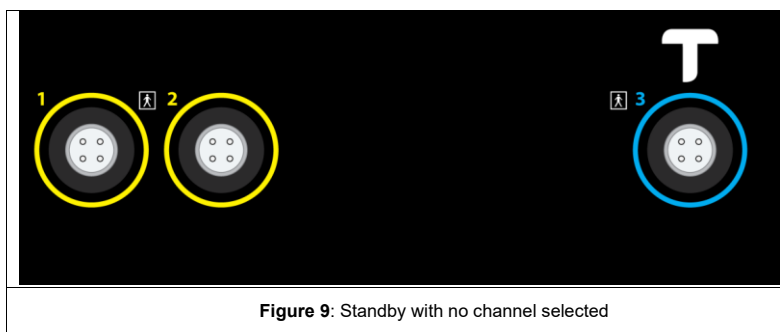
Power up the Generator

- Connect the mains lead to the rear panel of the ~~generator~~Generator.
- Press the Mains rocker switch up, "I", to switch on.
- The ~~generator~~Generator will now power up with an audio indicator flourish of tones, conduct a brief full screen illumination test as shown in ~~Figure 8~~Figure 8Figure 8Figure 8Figure 8.

Formatted: Font: +Body (Calibri), 11 pt



- After approx one second the display will show:



Optional - Attach the "Cement" footswitch to the Rear Panel

- The tube with the **Yellow** + "1" collar attaches to the similarly marked connection.
- The tube with the **Yellow** + "2" collar attaches to the similarly marked connection.

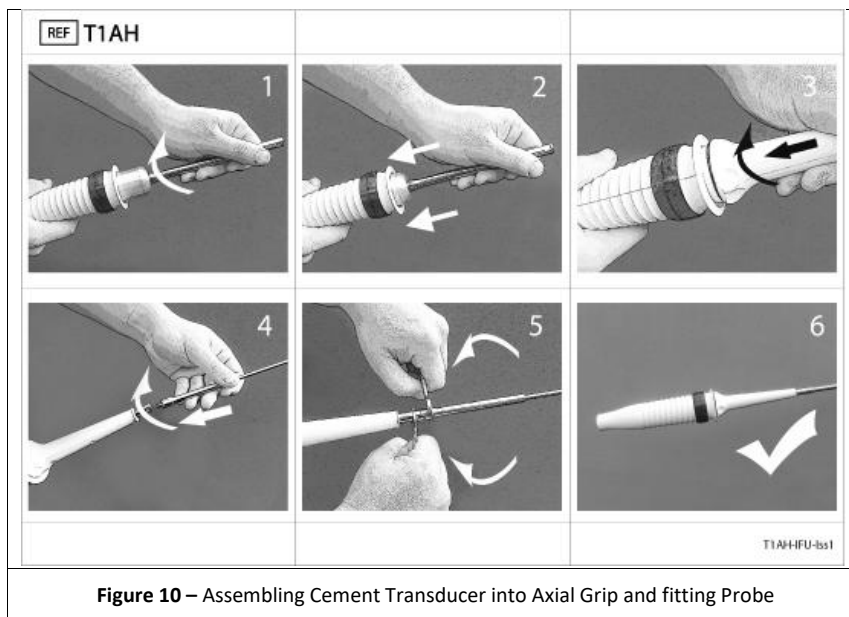
Optional - Attach the "Soft Tissue" footswitch to the Rear Panel.

- The tube with the **Grey** + "3" collar attaches to the similarly marked connection.
- The tube with the **Blue** + "3" collar attaches to the similarly marked connection.

◆ Alternatively, the finger switches may be used for the Soft Tissue Transducer.

Assemble the Cement Transducer + Probe (Applied Part)

- Select the pre-sterilized Cement **T**ransducer.
- Select the appropriate pre-sterilized reusable Probe.
- Select the Single Use Cement **H**andpiece (T1AH Axial Grip and Sleeve).
- Fit the Transducer into the Handpiece. Rotate and click into place – ensure that it is securely connected. See (1) & (2) in Figure 10 below.
- Slide the Handpiece **S**leeve over the Transducer horn. Locate the lugs at the base of the sleeve into the grip Handpiece and twist / click into place. See (3) below.
- Using supplied spanners, attach the Probe to the Transducer (4) & (5) below.



Optional - Attach the Soft Tissue Transducer (Applied Part)

- Remove the plug cap by pulling the knurled rings on the plug and the autoclave cap.
- Align the red dot on the metal plug with the red dot on the **generator** socket.
- Connect the Transducer to Generator Channel 3.
- Push the plug into the socket until it clicks.
- To remove, pull back on the knurled part of the plug - **Do not** pull the cable.

Cement Transducer - Connect

- Connect the Cement **c**Cable (Yellow coding) to Generator Channel 1 or 2 (Yellow rings) by aligning the red dot on the metal plug with the red dot on the **G**enerator socket. Push the plug into the socket until it clicks.
- Connect the other end of the Cement **c**Cable (Plastic connector) to the Cement **h**Handpiece. Align the arrow on the plastic connector with the slot in the rear connector on the Axial **g**Grip. Click into place.
- To remove, pull back on the knurled part of the plug - **Do not** pull the cable.

Cement Transducer - Initialise

- When the surgeon is ready to use the Transducer they must initialize it to enable power. This is done simply by pressing the black activate button or yellow **footswitchFootswitch (Ch1 or Ch2)**. The display shows:



Figure 11: Cement + Channel 1 Selected

Cement Transducer - Activate

- Press the black activate button on the **H**andpiece or the yellow **footswitchFootswitch**. This will energize the Transducer and Probe for surgery. An audible indicator signifies power delivery with a continuous low pitched tone.



Figure 12: Channel 1, Cement Transducer Active

- To stop output release the **footswitchFootswitch** or activate button.
- If the **footswitchFootswitch** or activate button is held down continuously for 20 seconds the audio indicator sound will change from continuous to pulsed.
- After a further 5 seconds the **generatorGenerator** will terminate power.

- Output power is cut
- Audible indicator stops
- Displays Warning symbol + Transducer **S**Symbol + Clock symbol.
- These will remain illuminated as long as the **Footswitch**Footswitch or activate button pressed.
- When it is released, normal operation resumes

Cement Removal – Irrigation

It is important to irrigate the bone canal during the cement removal process for two reasons: - It clears the canal of cement debris, including dust, and it also helps maintain safe operating temperatures within the canal. Such Irrigation systems are standard provision in orthopaedic cases.

It is recommended that irrigation in the form of pulsed lavage or manual irrigation is used after each activation of any Probe.

Do not irrigate *during* an activation as this may compromise the efficacy of the ultrasound transmission, translating to a slower cement removal process. The TORS probes are designed to produce rapid heating at the Probe / PMMA interface, thus softening the cement. Any cooling during this process will impede the desired effect.

Any waste products should be disposed of as per normal hospital procedures.

Fume Extraction

It is recommended that a suitable extraction system is used to remove smoke, aerosols and any associated odors.

Any waste products should be disposed of as per normal hospital procedures.

Soft Tissue Transducer - Connect

Plug a Soft **T**issue **T**ransducer into front Channel 3 (Blue ring) by aligning the red dots on the plug and connector receptacle. Push connector and click into place.

Soft Tissue Transducer - Initialize

- Press the Ch3 grey **Footswitch** to select the channel.
- Audio triple downwards tone.
- Note that it will default to **Low pPower** (single bar). The display shows:



Figure 13 Soft **T**issue **T**ransducer Initialize

Soft Tissue Transducer - Change Power Level

- If **High pPower** is required press the grey **Footswitch** again.
- Audio triple upwards tone. The display shows:



Figure 14 Soft **T**issue **T**ransducer - **high power Level-2** Selected

- Further presses will toggle between **Low** & **High pPower**.

Soft Tissue Transducer – Activate

- Press the Ch3 blue [footswitchFootswitch](#). This will energize the Soft [tTissue tT](#)ransducer.
- Audio continuous tone signifies power delivery. Display shows:

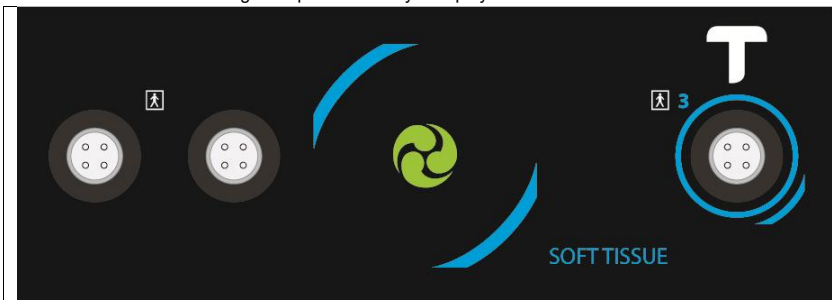
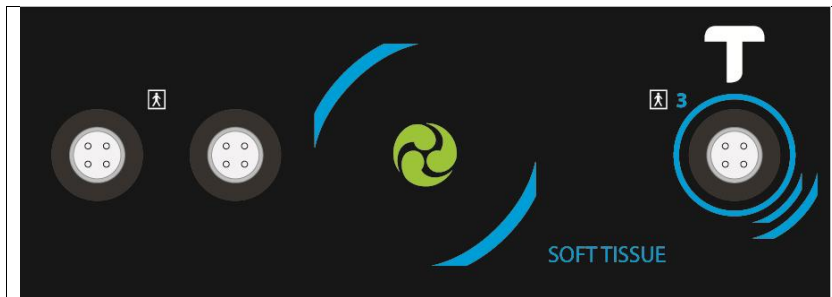


Figure 15 Soft [tTissue tT](#)ransducer [aActive](#) - [low pPower](#) Level-1

or



Soft [tTissue tT](#)ransducer [aActive](#) - [high pPower](#) Level-2

Depending upon power level selected





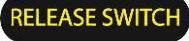

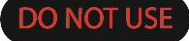


- To stop output release the [blue footswitchFootswitch](#) or [activate button](#).
- If the [footswitchFootswitch](#) or [activate button](#) is held down continuously for 20 seconds the audio indicator sound will change from continuous to pulsed.
- After a further 5 seconds the [generatorGenerator](#) will terminate power.
 - Output power is cut
 - Audible indicator stops
 - Displays Warning symbol + Transducer [sSymbol](#) + Clock symbol.
- These will remain illuminated as long as the [blue footswitchFootswitch](#) or [activate button](#) pressed.
- When it is released, normal operation resumes.


Turn Off TORS


- Press the black ON/OFF power switch at the rear of the unit. The screen will become entirely black.

Fault / Warning Indicators: (see also [Appendix 4: Tones and Beeps](#))


Formatted: Font: +Headings (Cambria), 14 pt

 <p>The Illuminated yellow triangle indicates that a fault has occurred.</p>	It will always be accompanied by:	
	i. An illuminated symbol indicating the part of the equipment where the fault has occurred	
		
	ii. Illuminated text instructing the user e.g. Reduce Pressure	<p>Or</p>  <p>Or</p>  <p>Or</p>  <p>Or</p> 
	If warning relates to time, the clock symbol also illuminates	
If warning is temperature related the temp symbol also illuminates		

	<p>Cement Transducer has been used 50 times. A high-low tone will sound. Prevents use of the Transducer.</p>
---	--


<p>+</p> 	<p>Or.....</p> <p>Soft issue transducer lifetime is nearly or over 100%. A high-low tone will sound. Prevents use of the Transducer.</p>
--	---

<p>FINAL SURGERY</p> <p>This will only occur during use</p>	<p>Cement transducer has been used 49 times. A high-low tone will sound. Allows use of Transducer for the duration of the current procedure.</p> <p>Or.....</p> <p>Soft issue transducer lifetime is nearly or over 100%. A high-low tone will sound. Allows use of Transducer for the duration of the current procedure.</p>
--	---

	<p>INCORRECT TRANSDUCER</p> <p>Transducer has been inserted in wrong socket. A high-low tone will sound. Prevents use of the Transducer.</p>
---	---


<p>REDUCE PRESSURE</p>	<p>Too much force is being applied to the TDCR during activation. Reduce applied force to improve performance. Excessive force applied during activation can impede cutting performance.</p>
-------------------------------	--

<p>RELEASE SWITCH</p>	<p>Handpiece button or footswitch held in. A high-low tone will sound.</p>
------------------------------	--




	<p>Transducer frequency too low. A high-low tone will sound. Try again after 5- 10s.</p>
---	--

<p><u>Alarm – This symbol is not currently used in TORS.</u></p>
--




- Formatted: Font: 3 pt
- Formatted: Centered
- Formatted: Indent: Left: 1.31 cm

	
---	--

Formatted: Font: (Default) Arial, 9 pt

	<p>The Ggenerator has detected a problem with the connections to the Transducer. The warning triangle light will illuminateflash, and a triple high-high-high tone will sound and repeat. Switch OFF then ON to clear fault. Re-energise the Transducer away from the patient – if fault recurs then change the Transducer and return for service.</p>
<p>+</p>  	

Formatted: Font: (Default) Arial, 9 pt, Bold

	<p>The Ggenerator has detected a problem. The warning triangle will illuminatelight will flash, and a triple high-high-high tone will sound and repeat. Switch OFF then ON to clear fault. Switch gGenerator back on – if fault recurs then return the system for service.</p>
<p>+</p>  	

Formatted: Font: (Default) Arial, 17 pt

DECONTAMINATION

The Cement **h**Handpiece is supplied pre-sterilized and is designed for single patient use.

After use the Cement **h**Handpiece should be disposed of as per normal theatre protocol, ideally in a suitably sized biohazard container.

The Cement **h**Handpiece is not designed for cleaning and sterilization after use. Post-operative cleaning and sterilization processes have not been validated by the manufacturer and any attempt at such process could cause the Cement **h**Handpiece to malfunction.

The Transducers and reusable Probes can be reused until their total activation time reaches 50 procedures or 250 minutes. The lifetime is checked at boot time and during use. After the lifetime limit is reached a warning message is displayed upon the **generatorGenerator's** screen.

Cement probes are detached from Cement **t**Transducers for decontamination

Soft **t**issue probe is NOT detached from Soft **t**issue **t**ransducer. It is washed and sterilized on the transducer for reuse.

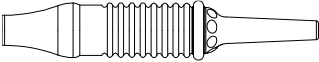




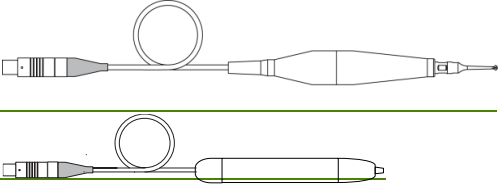





Formatted: Font: (Default) Arial, 9 pt

Formatted: Normal

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: Bold

Formatted: Space After: 0 pt

	<p>T1ASH - Cement aAxial gGrip hHandpiece (Single Use) Supplied Sterile (ethylene oxide) Do Not Use if packaging has been damaged!</p>	 
	<p>T1CT - Cement transducer (Reusable Applied Part)</p>	
	<p>T1TT - Soft tissue transducer + Probe (remains attached for decontamination and reuse) (Reusable Applied Part)</p>	
	<p>T1CC Cement rransducer cCable (Reusable)</p>	
	<p>T1FMU, T1SMU, T1PMU - Cement probes (Reusable Applied Part)</p>	

Formatted: Space After: 0 pt, Line spacing: single

Reprocessing Guidelines: Immediately After Use

- Immediately after use wipe down all components and remove any surplus body fluids and debris.

After Use

Detach the Handpiece from the Cement tTransducer by first removing the disposable sleeve. Then pull the Transducer and the HHandpiece apart.

Detach the Cement pProbe from the Cement tTransducer using supplied spanners.

Do not detached the Soft tTissue pProbe from the Soft tTissue tTransducer

Replace the tTransducers, HHandpiece cables, reusable pProbes and spanners back into their original wraps to reduce drying of soil while it is returned to the cleaning facility.

Single Use Items

Dispose of the Cement handpiece and sleeve as per hospital protocols.

Manual Cleaning

- Prepare enzymatic cleaning solution (e.g. Gigazyme Plus) to manufacturers instructions.
- Soak soiled instruments in enzymatic solution for 5 minutes.
- When cleaning, fully submerge the instruments in the cleaning solution. Brush with soft non-metallic bristle brush or cloth to remove all traces of blood and debris, concentrating on any crevices, seams, or other surface discontinuities. Clean holes and recesses using an appropriate brush ensuring that the full depth of the feature is reached. Ensure instruments are visibly clean before progressing to the next step.
- Rinse instruments thoroughly with clean running water for 2 minutes. Ensure that blind holes and recesses are repeatedly filled and emptied with running water.
- Dry instruments immediately after final rinse. Do not exceed 140 °C (285°F)

Note, currently no Automated Cleaning Process has been Validated and therefore only Manual Cleaning is authorized by the Manufacturer. Automated Cleaning

- ~~• Instruments may require manual cleaning prior to automated cleaning to improve the removal of adherent soil. Brush with a non-metallic bristle brush in enzymatic cleaning solution.~~
- ~~• Load instruments such that crevices, seams, surface discontinuities, holes and recesses can drain.~~
- ~~• Clean using the "Instruments" cycle in a validated washer-disinfector and a pH neutral cleaning agent intended for use in automatic cleaning (e.g. Getinge Enzymatic Detergent). The cleaning cycle should incorporate pre-rinse (minimum 3 minutes), wash (minimum 14 minutes), rinse (minimum 7 minutes), thermal-rinse (minimum 4 minutes at 93°C / 200°F, and drying steps — Do not exceed 140 °C (285°F).~~
- ~~• An alkaline cleaning solution such as, for example, Serchem pH Plus Detergent, with a pH up to 13.2 may be used instead of, or in addition to, an enzymatic solution.~~

Cleaning Inspection

- Inspect all instruments prior to sterilization or storage to ensure the complete removal of soiled surfaces.

- Visually inspect instruments, if soil is still present clean instruments again.
- Inspect cables for wear and damage, ensuring that no cracks, tears, or other damage is found.
- Check to see that Probes are free of scratches.
- Report any damage found to the Radley Scientific representative.

Wrapping

- Double wrap in accordance with local procedures, using standard procedures wrapping techniques such as those described in HTM 01-01 or ANSI/AAMI ST46-1993
- Label contents of wrapped tray using indelible marker or other sterilization compatible label system.

Sterilization

- Sterilization is best achieved on the day preceding the surgery, but must be at least one hour prior to use to allow the equipment to cool and stabilize.
- ⚠ Transducers are NOT to be submerged in water to expedite cooling.
- ⚠ Do NOT sterilize the ~~generator~~Generator or ~~footswitch~~Footswitches
- These components have been validated for sterilization by the following method in a vacuum autoclave. The parameters for this being 134-137°C (270°- 277°F) for a minimum of 3 and a maximum of 4 minutes. (If national standards dictate autoclaving up to 18 minutes, this is possible but not preferable.)

- ~~Next Use~~ - After decontamination, the Transducers, ~~Cement cables, probes and spanners~~, can be re-sterilized as above.

Formatted: Font: (Default) Arial, 10 pt

Formatted: Font: (Default) Arial, 10 pt

Formatted: Indent: Left: 0 cm, Hanging: 0.5 cm

End of Life Management

When a Transducer has reached the end of its lifetime please contact your local TORS supplier.

TORS generator and all reusable accessories (including ~~T~~transducers) should be recycled. Contact RSL for return instructions - see end page of this IFU.

A valid decontamination certificate must accompany each returned Transducer, ~~P~~probe or ~~Cement transducer~~ cable.

Opened in Error

In the event of a ~~Single Use disposable H~~handpiece being unwrapped in error, it may not be re-sterilized.

Care should be taken to ensure that the cables are not kinked during sterilization as this can produce cracks in the cable and reduce cable life.

Cleaning the Generator

The TORS Generator may be cleaned as follows:

- i. Dilute a neutral pH detergent according to the manufacturer's directions.
- ii. Using above solution, lightly moisten a soft, clean cloth. Wipe surfaces of the [generatorGenerator](#).
- iii. Using tap water, lightly moisten a soft, clean cloth. Wipe surfaces of the [generatorGenerator](#).
- iv. Dry [generatorGenerator](#) surfaces with a soft clean cloth.

Cleaning the Footswitch

The TORS Footswitches may be cleaned as follows:

- i. Dilute a neutral pH detergent according to the manufacturer's directions.
- ii. Using above solution, lightly moisten a soft, clean cloth. Wipe surfaces of the [feetswitchFootswitch](#).
- iii. Using tap water, lightly moisten a soft, clean cloth. Wipe surfaces of the [feetswitchFootswitch](#).
- iv. Dry [feetswitchFootswitch](#) surfaces with a soft clean cloth.
- v. Do not allow any water to enter into the air-hoses.
- vi. Do not detach the air hoses from the [feetswitchFootswitch](#).

Equipment Care

Probe Inspection

Before use it is recommended that the Probe be inspected for damage. TORS Probes are susceptible to damage if forced into contact with metal (e.g. hand instruments and clamps, etc.) when active. The consequence of scratching a Probe is to raise the mechanical stress in the region of the scratch (creating a "stress raiser") when the Probe is vibrating. If the stress is raised sufficiently there is a risk of the Probe suffering metal fatigue and cracking. Fatigue failure is more likely if the stress raiser is close to one of the fixed points of maximum mechanical stress in the vibrating Probe. If a Probe does suffer metal fatigue as a result of a stress raiser causing a crack, the Probe will no longer possess a resonant frequency "recognizable" by the generator. Report to the supplier any Transducer with a Probe that has a scratch, if the scratch is deep enough to be detected by sliding a fingernail over it.

Transducer / Probe Usage

The Cement and Soft tissue Transducers of TORS can be decontaminated and reused. This is an important factor in reducing the cost per case. However, they cannot be reused indefinitely and a lifetime is set within which optimum performance can be expected. In order to monitor the so-called "lifetime" of the either Transducer and Probe a memory chip is housed inside the cable plug of every Transducer.

For the Soft tissue transducer, every time the generator energizes the Transducer, the duration of the operation is monitored by the TORS then written to the memory chip, overwriting the previous total. When the total run time reaches pre-set markers for warning or termination, the TORS will alert the user.

For the Cement transducer, each time the Transducer is used for a procedure, the remaining number of allowed uses is monitored by the TORS then written to the memory chip, overwriting the previous total. When the total run reaches zero TORS will alert the user and further use will be prevented.

Storage of Equipment between Cases

It is recommended that the Transducers (with associated cables) are stored, between cases, in a large autoclave tray provided and again, care should be taken to ensure that the cables are not kinked close to the connector.

Maintenance

The hospital is responsible for ensuring that the unit has an electrical safety check performed by qualified service personnel at least once a year.

Do not remove the covers from TORS.

TORS generator does not require periodic calibration. If the generator detects an internal problem it will display a "Service Due" on the rear LCD. If this is seen, contact RSL to arrange repair.

There are no user-serviceable parts in TORS.

Any damage to the Transducers or cables should be reported and the components returned to the supplier at the earliest opportunity.




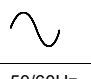
Technical Support

In line with the Company's policy of quality assurance and customer care, RSL Ltd. offers a telephone helpline service for users of TORS. This is provided by the TORS manufacturer, RSL Ltd at its facility in Devon, England.

HELPLINE Telephone No. as displayed on the unit is: +44 (0) 1364 653899

APPENDICES

Appendix 1: Markings on the TORS front and back panels

	Cement transducer selected
	Soft tissue transducer selected
	Cement transducer output - Channel 1
	Cement transducer output - Channel 2
	Soft tissue transducer output - Channel 3 (low power)
	Soft tissue transducer output - Channel 3 (high power)
	Transducer active
	Model n Number
	Serial n Number
	Manufacturer
	Date of m Manufacture
	FUSE rating - Time delay, 5A, High breaking capacity, 250 Volts AC, size 20mm
100-240V	Mains AC v Voltage range
	Alternating c Current
50/60Hz	Mains AC f Frequency
200VA	Input p Power
	Warning sign

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt






Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt




Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted Table

	Type BF e Equipment
36kHz	Cement f Frequency of o Output
<70W	Cement o Output p Power
For Intermittent Use OFF/ON/OFF 20/20s	Cement d Duty c Cycle
36kHz	Soft T issue f Frequency of o Output
<50W	Soft T issue o Output p Power
For Intermittent Use OFF/ON/OFF 3/30s	Soft T issue d Duty c Cycle
	Volume
	Footswitch connection
	Equipotential e Connection
	Output channel

Formatted Table

	<p>Follow instructions for use</p>
	<p>Electrical and Electronic equipment. Return waste to a collection system or treatment and recycling facilities. Follow decontamination instructions before returning waste.</p>
	<p>General Caution Sign</p>

Appendix 2: TORS Parts List

ITEM	Product Code
Cement Axial Grip Handpiece (Single Use)	T1AH
Cement 200 Ø9 Fan Probe (Reusable)	T1FM-9-200
Cement 300 Ø9 Fan Probe (Reusable)	T1FM-9-300
Cement 200 Ø15 Fan Probe (Reusable)	T1FM-15-200
Cement 300 Ø15 Fan Probe (Reusable)	T1FM-15-300
Cement 200 Ø20 Fan Probe (Reusable)	T1FM-20-200
Cement 300 Ø20 Fan Probe (Reusable)	T1FM-20-300
Cement 200 Ø8 Scraper Probe (Reusable)	T1SM-8.5-200
Cement 300 Ø8 Scraper Probe (Reusable)	T1SM-8.5-300
Cement 200 Ø6 Piercer Probe (Reusable)	T1PM-6.5-200
Cement 300 Ø6 Piercer Probe (Reusable)	T1PM-6.5-300
Cement 200 Ø8 Piercer Probe (Reusable)	T1PM-8.5-200
Cement 300 Ø8 Piercer Probe (Reusable)	T1PM-8.5-300
Cement 200 Ø8 Piercer/Scraper Probe (Reusable)	T1P/SM-8.5-200
Cement 300 Ø8 Piercer/Scraper Probe (Reusable)	T1P/SM-8.5-300
Cement 200 Ø8 Two-Lobe Groover Probe (Reusable)	T1GM-8.5-200
Cement 300 Ø8 Two-Lobe Groover Probe (Reusable)	T1GM-8.5-300
Cement 200 Ø9 Extractor Probe (Reusable)	T1EM-9.5-200
Cement 300 Ø9 Extractor Probe (Reusable)	T1EM-9.5-300
Soft Tissue Transducer + Probe (Reusable)	T1TT
Cement Transducer (Reusable)	T1CT
Cement Cable (Reusable)	T1CC
Generator	T1G
Generator Cart	T1C
Generator Carry Case	T1GC
Cement Footswitch	T1FC
Soft Tissue Footswitch	T1FT
Mains Cable	934
Autoclave Tray	T1AT
Autoclave Tray Mat	T1SM
Instrument Carry Case	T1IC

Appendix 3: Electromagnetic Interference

Table 1

Guidance and manufacturer's declaration – electromagnetic emissions		
TORS T1GLG4 is intended for use in the electromagnetic environment specified below. The customer or the user of TORS should assure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	TORS uses RF energy only for its internal function. Therefore its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class A	
Harmonic emissions IEC 61000-3-2	No testing – not connection to public mains network	TORS is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations /flicker emissions IEC 61000-3-3	No testing – not connection to public mains network	

Table 2

Guidance and manufacturer's declaration – electromagnetic immunity			
TORS is intended for use in the electromagnetic environment specified below. The customer or the user of TORS should assure that it is used in such an environment.			
IMMUNITY test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±6 kV contact ±8 kV air (1)	Floors should be conductive. No synthetic material should be used in the environment. The relative humidity should be in the range 40% to 60%.
Electrical fast transient/burst IEC61000-4-4	±2 kV for power supply lines	±2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode ± 2 kV common mode	±1 kV differential mode ± 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % U_T (2) (100 % dip in U_T) for 0.5 cycle at: 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315°	0 % U_T (100 % dip in U_T) for 0.5 cycle at: 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315°	Mains power quality should be that of a typical commercial or hospital environment. If the user of TORS requires continued operation during power mains interruptions, it is recommended that TORS be powered from an uninterruptible power supply.
	0 % U_T (100 % dip in U_T) for 1 cycles	0 % U_T (100 % dip in U_T) for 1 cycles	
	70 % U_T (30 % dip in U_T) For 25/30 cycles	70 % U_T (30 % dip in U_T) For 25/30 cycles	
	0 % U_T (100 % interrupt in U_T) for 250/300 cycles	0 % U_T (100 % interrupt in U_T) for 250/300 cycles	
Power frequency (50/60 Hz) magnetic field IEC61000-4-8	3 A/m	No testing	No magnetically sensitive components.
<p>NOTES: Mitigation applied because of environment. U_T is the a.c. mains voltage prior to application of the test level.</p>			

Table 3


Guidance and manufacturer's declaration – electromagnetic immunity			
TORS is intended for use in the electromagnetic environment specified below. The customer or the user of TORS should assure that it is used in such an environment.			
IMMUNITY test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3V rms 150kHz to 80 MHz Outside ISM bands	3V rms	Portable and mobile RF communications equipment should be used no closer to any part of TORS, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.2\sqrt{P}$ $d = 1.2\sqrt{P}$, 80MHz to 800MHz $d = 2.3\sqrt{P}$, 800MHz to 2.3GHz
Radiated RF IEC 61000-4-3	6V rms In ISM bands 0.15MHz to 80MHz 80% AM at 1kHz 3 V/m 80MHz to 2.5GHz	6V rms 3 V/m	where P is the maximum power output rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b
Immunity to proximity fields from RF wireless communications equipment	9 V/m 710MHz, 745MHz, 780MHz, 5240MHz, 5500MHz, 5785MHz 27 V/m 385MHz 28 V/m 450MHz, 810MHz, 870MHz, 930MHz, 1720MHz, 1845MHz, 1970MHz, 2450MHz	9V/m 27V/m 28V/m	Interference may occur in the vicinity of equipment marked with the following symbol. 
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
<i>a Field strengths from fixed transmitters , such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the in the location in which TORS is used exceed the applicable RF compliance level above, TORS should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating TORS.</i>			
<i>b Over the frequency range 150 kHz to 80 MHz , field strengths should be less than 3 V/m</i>			

Table 4

Recommended separation distances between portable and mobile RF communications equipment and TORS			
TORS is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of TORS can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment and TORS as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter	Separation distance according to frequency of transmitter		
	m		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800MHz to 2.5 GHz $d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance <i>d</i> in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where <i>P</i> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Appendix 4: Tones and Beeps

Rising flourish	TORS is booting to standby mode
No Tone	TORS is in standby mode
Triple tone – low-medium-low pitch	Handpiece a Activate b Button has been pressed to initialize handpiece (Cement) Handpiece t Toggle b Button has been pressed to initialize handpiece (Soft Tissue)
Continuous low pitched tone	Acoustic output at the LOW power level
Continuous high pitched tone	Acoustic output at the HIGH power level
Triple tone – low-medium-high pitch	Soft t issue t ransducer is changing from LOW to HIGH power
Triple tone – high-medium-low pitch	Soft t issue t ransducer is changing from HIGH to LOW power
Double beep –high>low pitch	1. The G enerator has reset itself after a minor problem such as a mistuning event, time-out or over-temperature warning or 2. The G enerator has encountered a minor problem that requires power switch to be cycled OFF/ON
Triple beep – high>high>high pitch	A more serious problem has been detected

Appendix 5: Rear Panel Messages

The LCD on the rear panel will display messages that indicate the status of the equipment. The following table shows the possible displays:

Message	Status	Action required
Radley Scientific Ltd TORS Issue x	Start up message. Shows the Issue number of the software "X"	Now awaiting Transducer to be plugged in
Cement / Soft Tissue Ready	The Toggle switch has been operated and the relevant Transducer is ready to be activated.	Now awaiting use
36500Hz Low Power	Low power is selected (STT Only) (via Soft tissue footswitch-or handpiece)	No action required
36500Hz High Power	High power is selected (STT Only) (via Soft tissue footswitch-or handpiece)	No action required
36000Hz Cement Active OR Soft Tissue Active	While Transducer is active the top line will show the frequency. After releasing the switch it will display the final running frequency.	No action required
Active Too long Release Switches	The active button has been held on for too long. No output.	Release activate button on handpiece or on footswitchFootswitch
Release Switches	Either a Activate or f Toggle button (or F ootswitches) has been pressed at switch on.	Release a Activate or f Toggle button when switching on Generator
Check Transducer Ease Grip and Retry	Transducer has been loaded too heavily Transducer is too hot.	Release switch, then reactivate using less pressure on the jaw Allow Transducer to cool.
Transducer Limit	Transducer is approaching 100% lifetime. Final use.	Replace Transducer after use
Change Tdcr Transducer Limit Transducer Limit Ready	Transducer has reached 100% lifetime and must be replaced.	Replace Transducer
Change Transducer Restart	Transducer frequency too low and feedback signal is low.	Switch supply off & on. Replace Transducer if seen 3 times
Transducer Leakage Change Tdcr	Generator has detected voltage on the Transducer.	Switch supply off & on. Replace Transducer if seen 3 times
Frequency Error Service due	The G enerator has detected a serious internal problem.	Switch supply off & on. If message seen again Generator requires service.

Appendix 6: Technical Specification

Model Nos:	See Appendix 2	
Dimensions:	Generator:	340mm (width) x 95mm (height) x 375mm (depth) 13.4" (width) x 3.7" (height) x 13.4" (depth)
Weight:	Generator:	7.6 kg
	Transport case:	13.8 kg (loaded with generator/Generator)
	Transducer:	0.37 kg
Fuse Type:	T5A, 250V, 20mm (2 off)	
Cordset	Contact RSL for recommended type	
Power supply input	100V - 240V, 50/60Hz	
Power consumption	200VA	
Cement Output - Frequency of operation	36kHz	
Cement Output – Power	<70W	
Cement Mode of Operation	Intermittent ON/OFF, 20/20s	
Cement Transducer Classification:	Type BF	
Soft Tissue Output - Frequency of operation	36kHz	
Soft Tissue Output – Power	<50W	
Soft Tissue Mode of Operation	Intermittent ON/OFF, 3/30s	
Soft Tissue Transducer Classification:	Type BF	
Accuracy of frequency display	1%	
Insulation Classification	Generator:	Class 1
Transducer	Titanium, s Stainless steel, and plastic Autoclavable maximum 50 cycles/	
Environment for Transportation & Storage:	Temperature: -10°C to +50°C Relative humidity: 10% to 90% Atmospheric pressure: 50kPa to 106kPa	
Environment for Use:	Temperature: +10°C to +30°C Relative humidity: 30% to 75% Atmospheric pressure: 81kPa to 106kPa	

RSL Ltd. will make available on request circuit diagrams, component part lists, descriptions, calibration instructions to assist to service personnel in parts repair.

TORS has been designed and built in accordance with ISO 13485: 2016 Quality Assurance standard for medical devices.

CE conformance has been certified and the equipment complies with:

IEC 60601-1:2005 + CORR. 1:2006 + CORR.2:2007 +
A1:2012
EN 60601-1:2006 + A11:2011 + A1:2013
ANSI/AAMI ES60601-1:2005/ (R) 2012

Appendix 7: Warranty Statement

Subject to the terms and conditions listed below, Radley Scientific Ltd. (hereafter called "the Company") guarantees to replace or repair free of charge any defective parts of TORS notified within the warranty period. This applies to the hardware defined below for the purposes of warranty claims made by any party supplied directly by Company.

Conditions of Warranty:

1. The Warranty covers the electronic system and ancillary components including the ~~generator~~Generator and ~~footswitch~~Footswitch, and applies to defects or faults that are reported to the Company, or its approved service facility, within ONE YEAR from the date of purchase.
2. Where transducers, cables or reusable probes are purchased separately or as part of the electronic system, the Warranty applies within the criteria of Clauses 5 and 6.
3. Where transducers, cables or reusable probes remain the property of the manufacturer the Warranty does not apply. If the Company loans a handset for use with the rest of the system the Company retains responsibility for performance.
4. The Warranty specifically excludes the disposables (~~H~~Handpieces and ~~S~~Single ~~U~~use probes). The disposables are detachable consumables for use with the ~~T~~Transducer. They are for ~~S~~Single ~~U~~use and are disposed of thereafter. Its condition and sterility is guaranteed up to the expiry date shown on the packaging label.
5. The Warranty relates to any components in Clauses 1 and 2, which are faulty or below specification as a result of inferior workmanship or materials. At the Company's expense the Company will conduct analysis, repair such faults, update to the relevant specification, re-test, re-calibrate and certify as necessary. For handsets this will take into account the additional criteria of Clause 6.
6. Transducer warranty criteria are defined as follows.
 - Cable, connector, acoustic system within (including waveguide) and transducer case are guaranteed against defective parts, materials or workmanship for 9 months from invoice date whether used or not, or until the transducer reaches end of life ("transducer limit") through use.
7. This Warranty does not cover accidental damage to, or fair wear and tear on, any part of the system.
8. This Warranty is only valid for use in the purchasing hospital and when used by personnel trained by the Company or its authorised distributors.
9. If a repair of a fault is required following expiry of the warranty, the Company will conduct analysis and submit a quotation to the customer stating recommended action and cost to the customer.
10. In the unlikely instance of the equipment requiring recall the Company will contact the purchasing hospital immediately and will, at the Company's expense, collect the relevant components.



Bremridge House,

Bremridge,
Ashburton

S. Devon TQ13 7JX
UK

Tel: (44) (1364) 653899 – Helpline

www.radleyscientific.co.uk

