

TORS

Torsional Orthopaedic Revision System

ARCHIVED COPY.

User Manual

DO NOT USE



CONTENTS

CONVENTIONS USED IN THIS DOCUMENT 755

Warning:..... 755

Caution:..... 755

CONTACT INFORMATION..... 755

INTENDED USE 866

PLEASE NOTE:..... 866

INDICATIONS..... 866

CONTRAINDICATIONS..... 866

THE TORS SYSTEM..... 977

 GENERATOR..... 977

Front Panel..... 977

Rear Panel..... 977

TORS Transducers & Cables..... 1088

SAFETY  124010

 WARNINGS..... 124010

 CAUTIONS 144111

 COMPLICATIONS AND POTENTIAL SIDE-EFFECTS 164212

INSTRUCTIONS FOR USE..... 184414

 SETTING UP THE TORS SYSTEM 184414

Power up the Generator..... 184414

Attach the "Cement" footswitch to the Rear Panel..... 184414

Optional - Attach the "Soft Tissue" footswitch to the Rear Panel..... 184414

Assemble the Cement Transducer + Probe (Applied Part)..... 194515

Optional - Attach the Soft Tissue Transducer (Applied Part)..... 194515

Cement Transducer - Connect..... 214616

Cement Transducer - Initialise..... 214616

Cement Transducer - Activate..... 214616

Cement Removal – Irrigation..... 224717

Fume Extraction..... 224717

ARCHIVED COPY:
DO NOT USE


<i>Soft Tissue Transducer - Connect</i>	241818
<i>Soft Tissue Transducer – Initialize</i>	241818
<i>Soft Tissue Transducer – Change Power Level</i>	241818
<i>Soft Tissue Transducer – Activate</i>	261919
<i>Turn Off TORS</i>	261919
FAULT / WARNING INDICATORS: (SEE ALSO APPENDIX 4: TONES AND BEEPS)	272020
DECONTAMINATION	322323
REPROCESSING GUIDELINES: IMMEDIATELY AFTER USE	342424
<i>After Use</i>	362624
<i>Single Use Items</i>	362624
MANUAL CLEANING	362624
AUTOMATED CLEANING	3424
AUTOMATED CLEANING	36262424
CLEANING INSPECTION	372624
WRAPPING	37262425
STERILIZATION	372625
<i>End of Life Management</i>	372625
<i>Opened in Error</i>	372625
<i>Cleaning the Generator</i>	392626
<i>Cleaning the Footswitch</i>	392626
EQUIPMENT CARE	402727
<i>Probe Inspection</i>	402727
<i>Transducer / Probe Usage</i>	402727
<i>Storage of Equipment between Cases</i>	402727
<i>Maintenance</i>	402727
<i>Technical Support</i>	402727
APPENDICES	412828
APPENDIX 1: MARKINGS ON THE TORS FRONT AND BACK PANELS	412828
APPENDIX 2: TORS PARTS LIST	4430313130
APPENDIX 3: ELECTROMAGNETIC INTERFERENCE	4631323231

ARCHIVED COPY:

DO NOT USE

Formatted: Font: (Default) Calibri, English (United States), Small caps

Formatted: Normal

APPENDIX 4: TONES AND BEEPS.....	4934353534
APPENDIX 5: REAR PANEL MESSAGES.....	5035363635
APPENDIX 6: TECHNICAL SPECIFICATION.....	5236373736
APPENDIX 7: WARRANTY STATEMENT.....	5337383837
CONVENTIONS USED IN THIS DOCUMENT.....	4
<i>Warning:</i>	4
<i>Caution:</i>	4
CONTACT INFORMATION.....	4
INTENDED USE.....	5
PLEASE NOTE:.....	5
INDICATIONS.....	5
CONTRAINDICATIONS.....	5
THE TORS SYSTEM.....	6
GENERATOR.....	6
<i>Front Panel</i>	6
<i>Rear Panel</i>	6
<i>TORS Transducers & Cables</i>	7
SAFETY 	9
WARNINGS.....	9
Cautions.....	10
INSTRUCTIONS FOR USE.....	11
SETTING UP THE TORS SYSTEM.....	11
<i>Power up the Generator</i>	11
<i>Attach the “Cement” footswitch to the Rear Panel</i>	11
<i>Optional— Attach the “Soft Tissue” footswitch to the Rear Panel</i>	11
<i>Assemble the Cement Transducer + Probe</i>	12
<i>Optional— Attach the Soft Tissue + Single Use Handpiece</i>	12
<i>Cement Transducer— Connect</i>	13
<i>Cement Transducer— Initialise</i>	13
<i>Cement Transducer— Activate</i>	13
<i>Cement Removal— Irrigation</i>	14

**ARCHIVED COPY:
DO NOT USE**

Formatted: Normal, Indent: First line: 0.42 cm

<i>Soft Tissue Transducer – Connect</i>	15
<i>Soft Tissue Transducer – Initialise</i>	15
<i>Soft Tissue Transducer – Change Power Level</i>	15
<i>Soft Tissue Transducer – Activate</i>	16
<i>Turn Off TORS</i>	16
FAULT / WARNING INDICATORS: (SEE ALSO APPENDIX 4: TONES AND BEEPS)	17
DECONTAMINATION	20
REPROCESSING GUIDELINES: IMMEDIATELY AFTER USE	21
<i>After Use</i>	21
MANUAL CLEANING	21
AUTOMATED CLEANING	21
CLEANING INSPECTION	21
WRAPPING	22
STERILIZATION	22
<i>Opened in Error</i>	22
<i>Cleaning the Generator</i>	23
<i>Cleaning the Footswitch</i>	23
EQUIPMENT CARE	24
<i>Probe Inspection</i>	24
<i>Transducer / Probe Usage</i>	24
<i>Storage of Equipment between Cases</i>	24
<i>Maintenance</i>	24
<i>Technical Support</i>	24
APPENDICES	25
APPENDIX 1: MARKINGS ON THE TORS FRONT AND BACK PANELS	25
APPENDIX 2: TORS PARTS LIST	27
APPENDIX 3: ELECTROMAGNETIC INTERFERENCE	28
APPENDIX 4: TONES AND BEEPS	31
APPENDIX 5: REAR PANEL MESSAGES	32
APPENDIX 6: TECHNICAL SPECIFICATION	33
APPENDIX 7: WARRANTY STATEMENT	34

**ARCHIVED COPY:
DO NOT USE**

|

**ARCHIVED COPY:
DO NOT USE**

Conventions used in this document

Warning:



A statement that if not strictly followed could result in injury or loss of life, or a statement that is otherwise required to be a warning.

Caution:

A statement that if not strictly followed could result in damage to the equipment

Contact Information

If the unit requires repair please contact us through your local supplier or directly:

RSL Ltd.
Bremridge House
Ashburton
S. Devon
TQ13 7JX
UK

Tel: +44 (0) 1364 653899

www.radleyscientific.co.uk

**ARCHIVED COPY:
DO NOT USE**

Intended use

TORS is an ultrasonic surgical device intended to be used for removal of polymethylmethacrylate (PMMA) bone cement in orthopaedic applications. It is also intended to be used for cutting and cauterising soft tissue during skin and muscle incision to access bone joints in limbs.

Please note:

- This document does not reference surgical techniques.
- The safety & effectiveness of any ultrasonic surgical equipment is mostly dependent upon the surgeon and nursing staff.

Indications

TORS is indicated for use in the removal of polymethylmethacrylate (PMMA) bone cement in orthopaedic applications. It is also indicated for soft tissue surgical incisions when bleeding control and minimal thermal injury are important.

Contraindications

- Do not use, if in the judgment of the surgeon, the use of ultrasonic surgical techniques are not in the best interest of the patient.
- Do not use for incising bone.
- The Soft Tissue Transducer is not indicated for the cutting of vessels exceeding 1.5mm in diameter or bone tissue.

**ARCHIVED COPY:
DO NOT USE**

Please read all Warning & Cautions contained in this document

The TORS System

Generator

Front Panel

The screen is highly visible, especially in a darkened operating theatre, indicates operating mode selected and gives concise instructions on action required in the event of disruption to normal operation.



Figure 1: TORS Generator Front Panel

Rear Panel

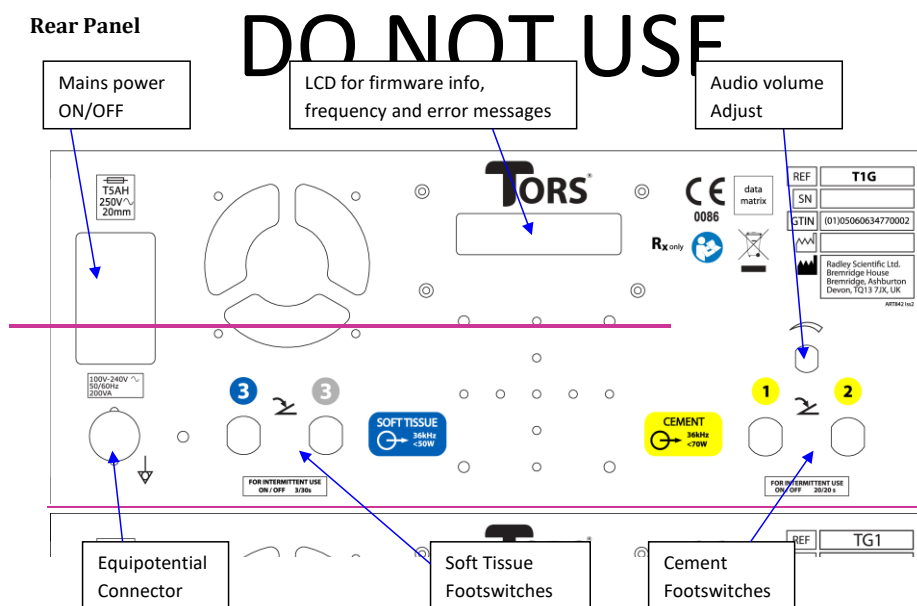


Figure 2: TORS Generator Rear Panel

There is only one user control on the generator itself:

- Volume control (on the rear panel)

TORS Transducers & Cables

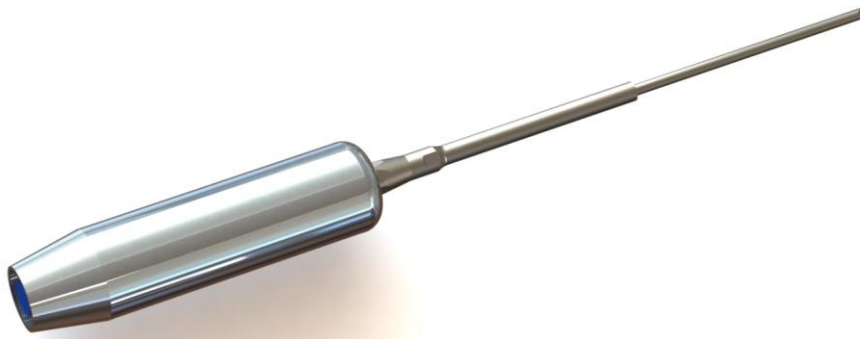


Figure 3: TORS Cement Transducer + Probe



Figure 4: TORS Soft Tissue Transducer + Probe

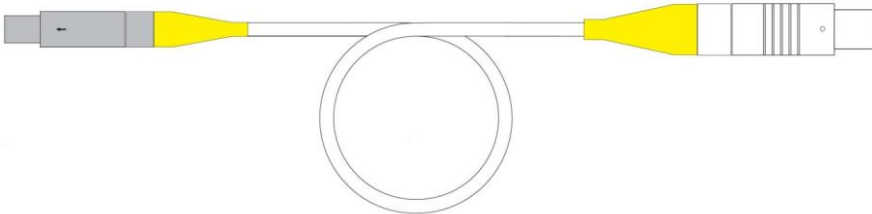


Figure 5: Cement Transducer Cable

|



Figure 6 – Axial Grip

ARCHIVED COPY:

DO NOT USE



Figure 7 – Axial Grip with TDCR and Probe fitted

Safety

Warnings

- This device should only be used by surgeons who are (1) trained in the types of surgical procedures that are to be carried out and (2) trained in the specific use of ultrasonic surgical instruments.
- Do not operate TORS in a potentially explosive or flammable area, or in oxygen rich environments. Note that if the Probe touches anything metal while activated, sparks may be produced.
- The use of TORS in cases where a patient or an operator of the unit has been fitted with a cardiac pacemaker is left to the discretion of the consultant in charge of the procedure although no electrical flow is made through or into the patient and many clinical papers on the use of ultrasound in the proximity of cardiac pacemakers have been published.
- Appropriate protective measures, including smoke evacuation, should be taken to protect users from any smoke, or other aerosols, produced by the use of surgical ultrasound.
- Use extreme caution when using in the vicinity of nerves.
- Handle the acoustics (Transducers and Probes) with care.
- Do not attempt to modify the acoustics. No modification of this equipment is allowed.
- The use of transducers or handpieces not supplied as part of the TORS system may damage the generator and create a safety hazard for the operator and patient.
- Avoid touching or holding the blade at the end of the energized Probe. It is designed to produce heat!
- Other than the activated tip do not allow probe to be in contact with tissue.
- Irrigate the venous canal after each activation.
- If the probe is deactivated inside heated PMMA it is possible the cooling cement will adhere very strongly to the probe, in which case it may become difficult to reactivate the probe in order to remove it from the cement.
- When guiding the TORS Cement Probe through PMMA, care must be taken not to exert excessive force in order to maintain control of the direction of travel of the Probe.
- Use of excessive force in vicinity of thin bone may cause perforation.
- Avoid resting the blade on skin or other tissue for at least 10s after cessation of energizing as it will have become hot while cutting. Ultrasonic systems dissipate heat quickly. (User may quench hot Probe tip / blade in saline)
- Care should be taken when in contact with tissue between activations, in case accidental activation should occur.
- If the Probe comes into contact with bone, using recommended cooling strategies a bone temperature rise of up to 5^oC may be seen.
- Allowing active tip to contact bone may cause damage to the bone.
- Wherever possible avoid contact between the side of active Probe and patient tissue.
- Mains isolation is achieved by use of the double pole switch located on the rear panel. DO NOT position the equipment -in such a way as to make access to this disconnection switch difficult.
- To avoid the risk of electric shock, this equipment must only be connected to a mains supply with protective earth.
- TORS is not MR safe or MR compatible.
- Do not re-sterilize or reuse any Single Use parts.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of TORS, including cables specified by RSL. Otherwise, degradation of the performance of this equipment could result.
- TORS should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary TORS should be observed to verify normal operation in the configuration in which it will be used
- The use of accessories, transducers and cables other than those specified may result in increased emissions, decreased immunity or improper operation.
- There are no user serviceable parts.

ARCHIVED COPY:
DO NOT USE

Cautions

- U.S. law restricts this device to sale by or on the order of a licensed physician.
- The nursing staff must be trained so that they are familiar with the equipment to be used.
- Avoid allowing an energized Probe to come into contact with any metal surface.
- Do not sterilize the TORS generator or ~~feetswitch~~Footswitch
- Do not block or otherwise restrict the vents on the rear and bottom panels.
- Never allow electrosurgical equipment to contact the Probe.
- Always transport TORS using the supplied transport cases.
- The TORS generator has an Equipotential Terminal on the back panel. This is provided for compatibility with other medical systems requiring such connection. This conductor is not intended for protective earthing. Refer to EN 60601-1 for details of use with the System.

- TORS:
 - i. should be used only for those procedures for which it is indicated.
 - ii. should be used with an appropriate power level commensurate with the required task.
 - iii. should be used with correct surgical technique.

If used correctly the TORS should cut and coagulate tissue and remove PMMA cement to the satisfaction of the user.

The equipment should inform the user of its status either by audio or visual means, unless ultrasound output is not possible. However, if a fault with the audio or visual indicators occurs whilst a cut is in progress, it may allow that cut to be completed without indication.

• Electromagnetic Interference

- This equipment is only suitable for use in hospital operating theatre / operating room.
- The performance of TORS may be degraded if it is subjected to electromagnetic disturbances e.g. an incorrect display that clears automatically once activate switch is pressed/depressed.
- Do not use TORS simultaneously with laser equipment or high frequency surgical equipment.
- This equipment has been tested and found to comply with the limits for a medical device.

However should interference occur, the user can try the following measures:-

1. Turn equipment off and on to confirm the source of the interference.
2. Increase separation between this equipment and other devices.
3. Connect this equipment to a power socket different from that to which the other devices are connected.
4. Consult medical physics department.

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Add space between paragraphs of the same style

ARCHIVED COPY.
DO NOT USE

Formatted: Font: (Default) Arial, 9 pt, Bold

Formatted: Font: (Default) Arial, 9 pt

Formatted: Normal, No bullets or numbering

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Add space between paragraphs of the same style

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Space After: 10 pt, Add space between paragraphs of the same style

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Add space between paragraphs of the same style

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Space After: 10 pt, Add space between paragraphs of the same style

Formatted: Indent: Left: 0.63 cm, Space After: 3 pt

Formatted: Indent: Left: 0.63 cm, Hanging: 0.63 cm, Tab stops: Not at 0.63 cm

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Space Before: 10 pt, After: 10 pt

- The EMISSIONS characteristics of this equipment make it suitable for use in industrial areas and hospitals (CISPR 11 class A). If it is used in a residential environment (for which CISPR 11 class B is normally required) this equipment might not offer adequate protection to radio-frequency communication services. The user might need to take mitigation measures, such as relocating or re-orienting the equipment.
- Where replacement Transducers and Cables are required, these must be manufactured by RSL (see Appendix 2 for full listing),
- See also [Appendix 3: Electromagnetic Interference](#)~~Appendix 3: Electromagnetic InterferenceAppendix 3: Electromagnetic InterferenceAppendix 3: Electromagnetic InterferenceAppendix 3: Electromagnetic InterferenceAppendix 3: Electromagnetic Interference~~.

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Add space between paragraphs of the same style

Formatted: Indent: Left: 0 cm, Hanging: 0.63 cm, Space After: 10 pt, Add space between paragraphs of the same style

Formatted: Font: (Default) Arial, 9 pt

Formatted: Indent: Left: 0.63 cm, No bullets or numbering

ARCHIVED COPY:
DO NOT USE

Complications and Potential Side-effects

General risks and complications of arthroplasty surgery

- Bone loss during prosthesis and/or cement removal
- Cortical perforation
- Fracture in the bone around the artificial joint
- Recovery time, pain and potential for arthritis associated with deliberate osteotomy
- Nerve/vascular lesions
- Rebleeding or haematoma / seroma
- Ligament, artery or nerve damage in the area around the knee joint
- Blood clots or deep vein thrombosis or pulmonary embolism
- Haemorrhage caused by tissue injury or ineffective haemostasis
- Bone/joint infection
- Late infection
- Infection of the wound and/or the tissue around the artificial joint
- Thromboembolism
- Implants or metal parts remaining in joint
- Fatigue fracture
- Excess bone forming around the artificial joint (knee) restricting movement
- Excess scar tissue forming and restricting movement
- Kneecap or hip dislocation
- Numbness in the vicinity of the wound scar
- Loosening of the artificial joint (normally after 10-15 years)
- Joint wear and tear
- Joint stiffening
- Bone cement implantation syndrome
- Hypoxia
- Hypotension
- Cardiac arrhythmias
- Cardiovascular collapse
- Thermal injury due to cement polymerization (setting) reaction
- Contact inhalational exposure to PMMA monomer leading to hypersensitivity
- Asthmatic reactions
- Neurological symptoms
- Localized irritation
- Personal risks from anaesthetic/anaesthetic or the surgery itself

Potential side-effects of ultrasonic arthroplasty surgery

- Necrosis and clinical dysfunction in femoral cortex and surrounding soft tissue due to heat transmission from ultrasonic probe/cement interface
- Intramedullary canal perforation

Formatted: Font: (Default) Arial, 10 pt

Formatted: Font: (Default) Arial, 10 pt

Formatted: Font: 10 pt

Formatted: Font: 10 pt

Formatted: Indent: Left: 0.63 cm, No bullets or numbering

Formatted: Line spacing: 1.5 lines

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Formatted: Line spacing: 1.5 lines

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Formatted: Line spacing: 1.5 lines

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

ARCHIVED COPY:
DO NOT USE

Formatted: Line spacing: 1.5 lines

Formatted: Indent: Left: 1.9 cm, Line spacing: 1.5 lines, No bullets or numbering

Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering

Formatted: Line spacing: 1.5 lines

Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering

Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm

Formatted: Font: 10 pt

Formatted: Font: 10 pt

- Different tissues and bone cement will be heated to different degrees with the same ultrasound parameters
- Each cell type has a different susceptibility to thermal injury
- Pathologic humeral fracture
- Microscopic bone cracks
- Plume production (95% water, 5% cell debris)
- Radial nerve palsy
- Risk of cross-infection from re-sterilized reusable parts of the system

General complications of powered dissection in open soft tissue surgery

- Acute or delayed bleeding (blood loss), organ perforation or fistula formation
- Haemorrhage caused by tissue injury or ineffective haemostasis
- Necrosis due to excessive heat transmission from the dissector
- Infections at the site of application
- Nerve lesions/damage

ARCHIVED COPY:

Potential side-effects of ultrasonic powered dissection in open soft tissue surgery

- Necrosis in soft tissues due to cavitation shock waves or excessive heat from ultrasonic probe
- Visibility reduction due to plume production (95% water, 5% cell debris) or smoke production
- Risk of transmission of mutagenic or toxic chemicals in the plume
- Risk of cross-infection from re-sterilized reusable parts of the system
- Device breakage and loss of components in the body cavity

DO NOT USE

- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Font: 10 pt
- Formatted: Indent: Left: 1.9 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Line spacing: 1.5 lines
- Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Line spacing: 1.5 lines
- Formatted: Line spacing: 1.5 lines, Bulleted + Level: 1 + Aligned at: 1.27 cm + Indent at: 1.9 cm
- Formatted: Font: (Default) Arial, 9 pt
- Formatted: Normal, Indent: Left: 1.27 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Indent: Left: 0.63 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Indent: Left: 0.75 cm, Line spacing: 1.5 lines, No bullets or numbering
- Formatted: Normal, Left, Indent: Left: 0.75 cm, Line spacing: 1.5 lines
- Formatted: Font: 10 pt
- Formatted: Line spacing: 1.5 lines
- Formatted: Font: 9 pt
- Formatted: Font: 9 pt
- Formatted: Font: (Default) Arial, 9 pt

◆ Alternatively, the finger switches may be used for the Soft Tissue Transducer.

Assemble the Cement Transducer + Probe (Applied Part)

- Select the pre-sterilized Cement **T**ransducer.
- Select the appropriate pre-sterilized reusable Probe.
- Select the Single Use Cement **h**andpiece (T1AH Axial Grip and Sleeve).
- Fit the Transducer into the Handpiece. Rotate and click into place – ensure that it is securely connected. See (1) & (2) in Figure 10 below.
- Slide the Handpiece **s**leeve over the Transducer horn. Locate the lugs at the base of the sleeve into the grip Handpiece and twist / click into place. See (3) below.
- Using supplied spanners, -attach the Probe to the Transducer (4) & (5) below.

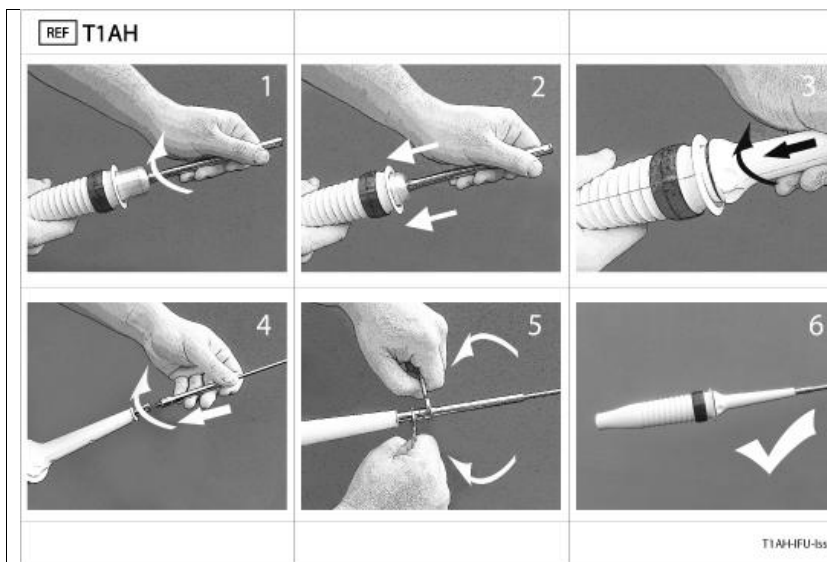


Figure 10 – Assembling Cement Transducer into Axial Grip and fitting Probe

Optional - Attach the Soft Tissue Transducer (Applied Part)

- Remove the plug cap by pulling the knurled rings on the plug and the autoclave cap.
- Align the red dot on the metal plug with the red dot on the **generator** socket.
- Connect the Transducer to Generator Channel 3.
- Push the plug into the socket until it clicks.
- To remove, pull back on the knurled part of the plug - **Do not** pull the cable.

**ARCHIVED COPY:
DO NOT USE**

Cement Transducer - Connect

- Connect the Cement **c**Cable (-Yellow coding-) to Generator Channel 1 or 2 (-Yellow rings-) by aligning the red dot on the metal plug with the red dot on the **G**enerator socket. Push the plug into the socket until it clicks.
- Connect the other end of the Cement **c**Cable (-Plastic connector-) to the Cement **h**Handpiece. Align the arrow on the plastic connector with the slot in the rear connector on the Axial **g**rip. Click into place.
- To remove, pull back on the knurled part of the plug - **Do not** pull the cable.

Cement Transducer - Initialise

- When the surgeon is ready to use the **T**ransducer, they must initialize it to enable power. This is done simply by pressing the black activate button or yellow **f**ootswitch (**Ch1 or Ch2**). The display shows:

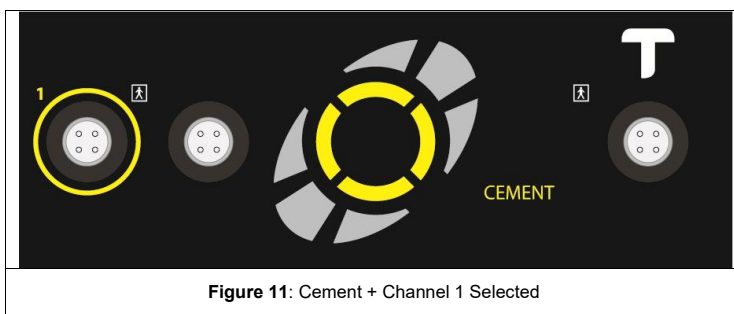


Figure 11: Cement + Channel 1 Selected

DO NOT USE

Cement Transducer - Activate

- Press the black activate button on the **H**andpiece or the yellow **f**ootswitch. This will energize the Transducer and Probe for surgery. An audible indicator signifies power delivery with a continuous low-pitched tone.

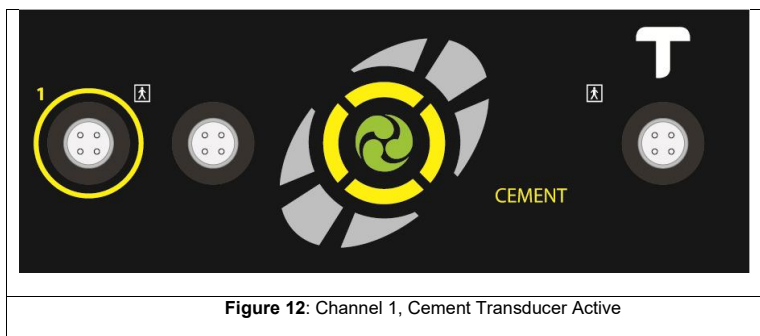


Figure 12: Channel 1, Cement Transducer Active

- To stop output, release the **f**ootswitch or activate button.
- If the **f**ootswitch or activate button is held down continuously for 20 seconds the audio indicator sound will change from continuous to pulsed.
- After a further 5 seconds the **g**enerator will terminate power.
 - Output power is cut
 - Audible indicator stops
 - Displays Warning symbol + Transducer **s**Symbol + Clock symbol.

- These will remain illuminated as long as the ~~footswitch~~Footswitch or activate button pressed.
- When it is released, normal operation resumes

Cement Removal - Irrigation

It is important to irrigate the bone canal during the cement removal process for two reasons: - It clears the canal of cement debris, including dust, and it also helps maintain safe operating temperatures within the canal. Such Irrigation systems are standard provision in orthopaedic cases.

It is recommended that irrigation in the form of pulsed lavage or manual irrigation is used after each activation of any Probe.

Do not irrigate *during* an activation as this may compromise the efficacy of the ultrasound transmission, translating to a slower cement removal process. The TORS probes are designed to produce rapid heating at the Probe / PMMA interface, thus softening the cement. Any cooling during this process will impede the desired effect.

Any waste products should be disposed of as per normal hospital procedures.

Fume Extraction

It is recommended that a suitable extraction system is used to remove smoke, aerosols and any associated odors.

Any waste products should be disposed of as per normal hospital procedures.

**ARCHIVED COPY:
DO NOT USE**

**ARCHIVED COPY:
DO NOT USE**

Soft Tissue Transducer - Connect

- Plug a Soft Tissue Transducer into front Channel 3 (-Blue ring-) by aligning the red dots on the plug and connector receptacle. Push connector and click into place.

Soft Tissue Transducer - Initialize

- Press the Ch3 grey footswitch to select the channel.
- Audio triple downwards tone.
- Note that it will default to Low Power (single bar). The display shows:

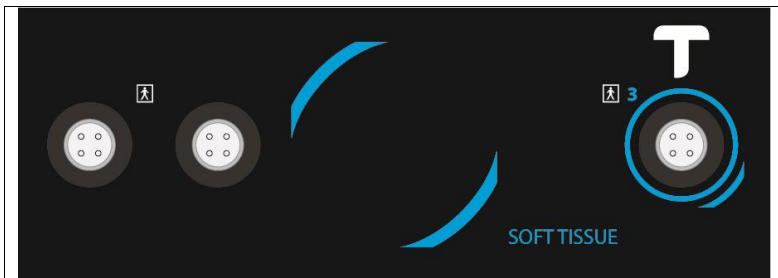


Figure 13 Soft Tissue Transducer initialize

ARCHIVED COPY:

Soft Tissue Transducer - Change Power Level

- If High Power is required press the grey footswitch again.
- Audio triple upwards tone. The display shows:

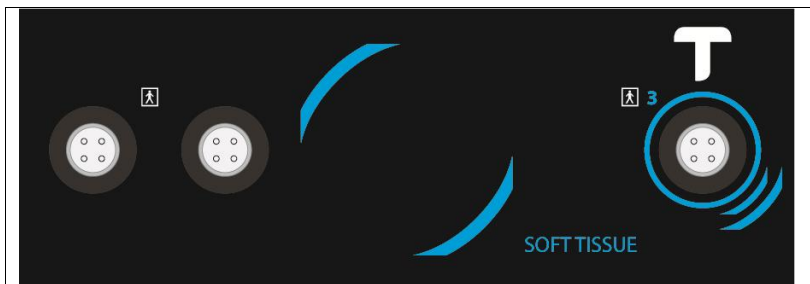


Figure 14 Soft Tissue Transducer - high power Level 2 Selected

- Further presses will toggle between Low & High Power.

|

**ARCHIVED COPY:
DO NOT USE**

Soft Tissue Transducer – Activate

- Press the Ch3 blue footswitchFootswitch. This will energize the Soft tTissue tTransducer.
- Audio continuous tone signifies power delivery. Display shows:


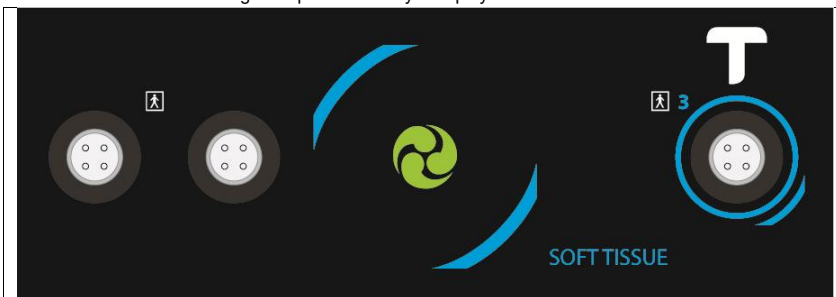


Figure 15 Soft tTissue ttransducer- aActive - -low pPower- Level-1

or

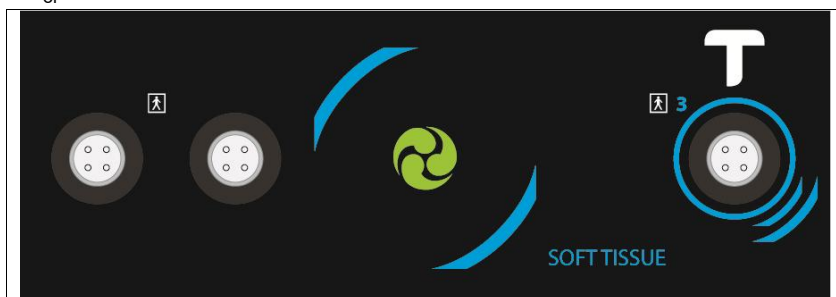

Soft tTissue ttransducer- aActive - -high pPower- Level-2

Depending upon power level selected










- To stop output, release the blue footswitchFootswitch-or-activate-button.
- If the footswitchFootswitch-or-activate-button is held down continuously for 20 seconds the audio indicator sound will change from continuous to pulsed.
- After a further 5 seconds the generatorGenerator will terminate power.
 - Output power is cut
 - Audible indicator stops
 - Displays Warning symbol + Transducer sSymbol + Clock symbol-
- These will remain illuminated as long as the blue footswitchFootswitch-or-activate-button pressed- .
- When it is released, normal operation resumes.


Turn Off TORS


- Press the black ON/OFF power switch at the rear of the unit. The screen will become entirely black.

Fault / Warning Indicators: (see also [Appendix 4: Tones and Beeps](#)~~Appendix 4: Tones and Beeps~~~~Appendix 4: Tones and Beeps~~~~Appendix 4: Tones and Beeps~~~~Appendix 4: Tones and Beeps~~)

Formatted: Font: +Headings (Cambria), 14 pt

 <p>The Illuminated yellow triangle indicates that a fault has occurred.</p>	It will always be accompanied by:	
	i. An illuminated symbol indicating the part of the equipment where the fault has occurred	
		
	ii. Illuminated text instructing the user e.g. Reduce Pressure	<p>Or</p>  <p>Or</p>  <p>Or</p>  <p>Or</p> 
	If warning relates to time, the clock symbol also illuminates	
If warning is temperature related the temp symbol also illuminates		

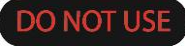	<p>Cement Transducer has been used 50 times. A high-low tone will sound. Prevents use of the Transducer.</p>
---	--


+	 <p>Or.....</p> <p>Soft tissue Transducer lifetime is nearly or over 100%. A high-low tone will sound. Prevents use of the Transducer.</p>
---	--

<p>FINAL SURGERY</p> <p>This will only occur during use</p>	<p>Cement Transducer has been used 49 times. A high-low tone will sound. Allows use of Transducer for the duration of the current procedure.</p> <hr/> <p>Or.....</p> <p>Soft tissue Transducer lifetime is nearly or over 100%. A high-low tone will sound. Allows use of Transducer for the duration of the current procedure.</p>
--	---

ARCHIVED COPY:


[Note. The working lifetime of the transducers and probes are detailed in the Decontamination section \(Page 23.\)](#)

Formatted: Indent: Left: 0 cm

	<p style="font-size: 24pt; font-weight: bold; text-align: center;">DO NOT USE</p> <p>INCORRECT TRANSDUCER Transducer has been inserted in wrong socket. A high-low tone will sound. Prevents use of the Transducer.</p>
---	--

<p>REDUCE PRESSURE</p>	<p>Too much force is being applied to the TDCR during activation. Reduce applied force to improve performance. Excessive force applied during activation can impede cutting performance.</p>
-------------------------------	--

<p>RELEASE SWITCH</p>	<p>Handpiece button or footswitch held in. A high-low tone will sound.</p>
------------------------------	---

	<p>Transducer frequency too low. A high-low tone will sound. Try again after 5- 10s.</p>
---	--

Formatted: Font: 3 pt

Formatted: Centered



Alarm – This symbol is not currently used in TORS.

Formatted: Font: (Default) Arial, 9 pt

Formatted: Indent: Left: 1.31 cm

ARCHIVED COPY:
DO NOT USE



The **G**enerator has detected a problem with the connections to the Transducer. The warning **triangle light** will **illuminate/flash**, and a triple high-high-high tone will sound and repeat. Switch OFF then ON to clear fault. Re-energise the Transducer away from the patient – if fault recurs then change the Transducer and return for service.

+



ARCHIVED COPY:
DO NOT USE

Formatted: Font: (Default) Arial, 9 pt, Bold



The **G**enerator has detected a problem. The warning **triangle will illuminate/light will flash**, and a triple high-high-high tone will sound and repeat. Switch OFF then ON to clear fault. Switch **g**enerator back on – if fault recurs then return the system for service.

+



Formatted: Font: (Default) Arial, 17 pt

|

**ARCHIVED COPY:
DO NOT USE**

DECONTAMINATION

The Cement **h**Handpiece is supplied pre-sterilized and is designed for single patient use.

After use the Cement **h**Handpiece should be disposed of as per normal theatre protocol, ideally in a suitably sized biohazard container.

The Cement **h**Handpiece is not designed for cleaning and sterilization after use. Post-operative cleaning and sterilization processes have not been validated by the manufacturer and any attempt at such process could cause the Cement **h**Handpiece to malfunction.

Both the Cement and Soft-Tissue Transducers can be reused until their total activation time reaches 25 procedures or 250 minutes. The lifetime is checked at boot time and during use. After the lifetime limit is reached a warning message is displayed upon the Generator's screen.






Reusable Cement Probes fall into two categories:

- T1FM-9-200, T1SM-8-200, T1CM-3-200, T1PM-6, 8, 10, 12, 14, 16, 18, 20-200 are all 5 minutes ON-time per case, over 50 cases = 250 minutes ON-time service life.
- T1SM-8, 10, 12, 16-200 all 10 minutes ON-time per case, over 50 cases = 500 minutes ON-time service life.

The Transducers and reusable Probes can be reused until their total activation time reaches 50 procedures or 250 minutes. The lifetime is checked at boot time and during use. After the lifetime limit is reached a warning message is displayed upon the generator's screen.

Cement probes are detached from Cement **t**ransducers for decontamination.

Soft **t**issue probe is NOT detached from Soft **t**issue **t**ransducer. It is washed and sterilized on the transducer for reuse.

	T1ASH - Cement aAxial-gGrip hHandpiece (Single Use) Supplied Sterile (ethylene oxide) Do Not Use if packaging has been damaged!	 
	T1CT - Cement tTransducer (Reusable Applied Part)	
	T1TT - Soft t issue tTransducer + P erobe (remains attached)	

Formatted: Space Before: 0 pt, After: 18 pt, Add space between paragraphs of the same style

Formatted: Heading 1, Space After: 18 pt

Formatted: Justified, Space After: 6 pt

Formatted: Justified, Space Before: 0 pt

Formatted: Justified, Space Before: 0 pt, After: 6 pt

Formatted: Justified, Space After: 6 pt, Add space between paragraphs of the same style

Formatted: Normal, Justified, Space After: 6 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt








Formatted: Font: (Default) Arial, 9 pt

Formatted: Justified, Space After: 6 pt, Line spacing: single

Formatted Table

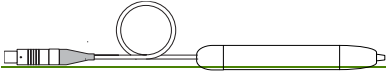Formatted: Space After: 0 pt, Line spacing: single

ARCHIVED COPY:
DO NOT USE

	<p><u>for decontamination and reuse</u> (Reusable Applied Part)</p>	
		
	<p>T1CC Cement Transducer Cable (Reusable)</p>	
	<p>T1FMU, T1SMU, T1PMU - Cement Probes (Reusable Applied Part)</p>	

Formatted: Centered

**ARCHIVED COPY:
DO NOT USE**

Reprocessing Guidelines: Immediately After Use

- Immediately after use wipe down all components and remove any surplus body fluids and debris.
- Activate the probe heads in the media within the cleaning cell, this will remove cement adhered to the surfaces.

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Font: 10 pt

After Use

Detach the Handpiece from the Cement transducer by first removing the disposable sleeve. Then pull the Transducer and the Handpiece apart.

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

Detach the Cement probe from the Cement transducer using supplied spanners.

Formatted: Indent: Left: 0 cm

Do not detached the Soft tissue probe from the Soft tissue transducer

Replace the Transducers, Handpiece cables, reusable Probes and spanners back into their original wraps to reduce drying of soil while it is returned to the cleaning facility.

Single Use Items

Dispose of the Cement handpiece and sleeve as per hospital protocols.

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Indent: Left: 0 cm

Reusable Items

Manual Cleaning

- Prepare enzymatic cleaning solution (e.g. Gigazyme Plus) to manufacturer's instructions.
- Soak soiled instruments in enzymatic solution for 5 minutes.
- When cleaning, fully submerge the instruments in the cleaning solution. Brush with soft non-metallic bristle brush or cloth to remove all traces of blood and debris, concentrating on any crevices, seams, or other surface discontinuities. Clean holes and recesses using an appropriate brush ensuring that the full depth of the feature is reached. Ensure instruments are visibly clean before progressing to the next step.
- Rinse instruments thoroughly with clean running water for 2 minutes. Ensure that blind holes and recesses are repeatedly filled and emptied with running water.
- Dry instruments immediately after final rinse. Do not exceed 140°C (285°F)

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Space Before: 12 pt, After: 6 pt

Automated Cleaning

- Instruments may require manual cleaning prior to automated cleaning to improve the removal of adherent soil. Brush with a non-metallic bristle brush in enzymatic cleaning solution.
- Load instruments such that crevices, seams, surface discontinuities, holes and recesses can drain.
- Clean using the "Instruments" cycle in a validated washer disinfectant and a pH neutral cleaning agent intended for use in automatic cleaning (e.g. Getinge Enzymatic Detergent). The cleaning cycle should incorporate pre-rinse (minimum 3 minutes), wash (minimum 14 minutes), rinse (minimum 7 minutes), thermal rinse (minimum 4 minutes at 93°C / 200°F, and drying steps - Do not exceed 140°C (285°F).
- An alkaline cleaning solution such as, for example, Serchem pH Plus Detergent, with a pH up to 13.2 may be used instead of, or in addition to, an enzymatic solution.

Formatted: Space Before: 12 pt, After: 6 pt

Cleaning Inspection

- Inspect all instruments prior to sterilization or storage to ensure the complete removal of soiled surfaces.
- Visually inspect instruments, if soil is still present clean instruments again.

ARCHIVED COPY:
DO NOT USE

- [Inspect cables for wear and damage, ensuring that no cracks, tears, or other damage is found.](#)
- [Check to see that Probes are free of scratches.](#)
- [Report any damage found to the Radley Scientific representative.](#)

Wrapping

- [Double wrap in accordance with local procedures, using standard procedures wrapping techniques such as those described in HTM 01-01 or ANSI/AAMI ST46-1993](#)
- [Label contents of wrapped tray using indelible marker or other sterilization compatible label system.](#)

Sterilization

- [Sterilization is best achieved on the day preceding the surgery, but must be at least one hour prior to use to allow the equipment to cool and stabilize.](#)
- ⚠ [Transducers are NOT to be submerged in water to expedite cooling.](#)
- ⚠ [Do NOT sterilize the Generator or Footswitches](#)
- [These components have been validated for sterilization by the following method in a vacuum autoclave. The parameters for this being 134-137°C \(270°- 277°F\) for a minimum of 3 and a maximum of 4 minutes. \(If national standards dictate autoclaving up to 18 minutes, this is possible but not preferable\)](#)
- [Next Use - After decontamination, the Transducers, Cement cables, probes and spanners can be re-sterilized as above.](#)

End of Life Management

[When a Transducer has reached the end of its lifetime please contact your local TORS supplier. TORS generator and all reusable accessories \(including Transducers\) should be recycled. Contact RSL for return instructions - see end page of this IFU.](#)

[A valid decontamination certificate must accompany each returned Transducer, Probe or Cement cable.](#)

Opened in Error

[In the event of a Single Use Handpiece being unwrapped in error, it may not be re-sterilized. Care should be taken to ensure that the cables are not kinked during sterilization as this can produce cracks in the cable and reduce cable life.](#)

Cleaning the Generator

[The TORS Generator may be cleaned as follows:](#)

- [Dilute a neutral pH detergent according to the manufacturer's directions.](#)
- [Using above solution, lightly moisten a soft, clean cloth. Wipe surfaces of the Generator.](#)
- [Using tap water, lightly moisten a soft, clean cloth. Wipe surfaces of the Generator.](#)
- [Dry Generator surfaces with a soft clean cloth.](#)

Cleaning the Footswitch

[The TORS Footswitches may be cleaned as follows:](#)

- [Dilute a neutral pH detergent according to the manufacturer's directions.](#)
- [Using above solution, lightly moisten a soft, clean cloth. Wipe surfaces of the Footswitch.](#)

ARCHIVED COPY:
DO NOT USE

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Indent: Left: 0 cm

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

Formatted: Font color: Auto

Formatted: Space Before: 12 pt, After: 6 pt

- iii. Using tap water, lightly moisten a soft, clean cloth. Wipe surfaces of the Footswitch.
- iv. Dry Footswitch surfaces with a soft clean cloth.
- v. Do not allow any water to enter into the air-hoses.
- vi. Do not detach the air hoses from the Footswitch.

Reprocessing Guidelines: Immediately After Use

- ~~Immediately after use wipe down all components and remove any surplus body fluids and debris.~~

Formatted: No bullets or numbering

After Use

~~Detach the Handpiece from the Cement TTransducer by first removing the disposable sleeve. Then pull the Transducer and the Hhandpiece apart.~~

~~Detach the Cement pProbe from the Cement TTransducer using supplied spanners.~~

~~Do not detached the Soft tTissue pProbe from the Soft tTissue tTransducer~~

~~Replace the Ttransducers, Hhandpiece cables, reusable Pprobes and spanners back into their original wraps to reduce drying of soil while it is returned to the cleaning facility.~~

Single Use Items

ARCHIVED COPY:

~~Dispose of the Cement handpiece and sleeve as per hospital protocols.~~

Manual Cleaning

DO NOT USE

- ~~Prepare enzymatic cleaning solution (e.g. Cegzyme Plus) to manufacturers instructions.~~
- ~~Soak soiled instruments in enzymatic solution for 5 minutes.~~
- ~~When cleaning, fully submerge the instruments in the cleaning solution. Brush with soft non-metallic bristle brush or cloth to remove all traces of blood and debris, concentrating on any crevices, seams, or other surface discontinuities. Clean holes and recesses using an appropriate brush ensuring that the full depth of the feature is reached. Ensure instruments are visibly clean before progressing to the next step.~~
- ~~Rinse instruments thoroughly with clean running water for 2 minutes. Ensure that blind holes and recesses are repeatedly filled and emptied with running water.~~
- ~~Dry instruments immediately after final rinse. Do not exceed 140 °C (285°F)~~

~~Note, currently no Automated Cleaning Process has been Validated and therefore only Manual Cleaning is authorized by the Manufacturer. Automated Cleaning~~

- ~~Instruments may require manual cleaning prior to automated cleaning to improve the removal of adherent soil. Brush with a non-metallic bristle brush in enzymatic cleaning solution.~~
- ~~Load instruments such that crevices, seams, surface discontinuities, holes and recesses can drain.~~
- ~~Clean using the "Instruments" cycle in a validated washer disinfecter and a pH neutral cleaning agent intended for use in automatic cleaning (e.g. Getinge Enzymatic Detergent). The cleaning cycle should incorporate pre-rinse (minimum 3 minutes), wash (minimum 14 minutes), rinse (minimum 7 minutes), thermal rinse (minimum 4 minutes at 93°C / 200°F, and drying steps – Do not exceed 140 °C (285°F).~~

- An alkaline cleaning solution such as, for example, Serchem pH Plus Detergent, with a pH up to 13.2 may be used instead of, or in addition to, an enzymatic solution.

Cleaning Inspection

- Inspect all instruments prior to sterilization or storage to ensure the complete removal of soiled surfaces.
- Visually inspect instruments, if soil is still present clean instruments again.
- Inspect cables for wear and damage, ensuring that no cracks, tears, or other damage is found.
- Check to see that Probes are free of scratches.
- Report any damage found to the Radley Scientific representative.

Wrapping

- Double wrap in accordance with local procedures, using standard procedures wrapping techniques such as those described in HTM 01-01 or ANSI/AAMI ST46-1993
- Label contents of wrapped tray using indelible marker or other sterilization compatible label system.

Sterilization

- Sterilization is best achieved on the day preceding the surgery, but must be at least one hour prior to use to allow the equipment to cool and stabilize.
- ⚠ Transducers are NOT to be submerged in water to expedite cooling.
- ⚠ Do NOT sterilize the generator Generator or footswitch Footswitches
- These components have been validated for sterilization by the following method in a vacuum autoclave. The parameters for this being 134-137°C (270°-277°F) for a minimum of 3 and a maximum of 4 minutes. (If national standards dictate autoclaving up to 18 minutes, this is possible but not preferable.)
- Next Use After decontamination, the Transducers, Cement cables, probes and spanners can be re-sterilized as above.

End of Life Management

When a Transducer has reached the end of its lifetime please contact your local TORS supplier.

TORS generator and all reusable accessories (including Ttransducers) should be recycled. Contact RSL for return instructions – see end page of this IFU.

A valid decontamination certificate must accompany each returned Transducer, Pprobe or Cement transducer cable.

Opened in Error

ARCHIVED COPY:
DO NOT USE

Formatted: Font: (Default) Arial, 10 pt

Formatted: Indent: Left: 0 cm, Hanging: 0.5 cm

In the event of a Single Use disposable Hhandpiece being unwrapped in error, it may not be re-sterilized.

Care should be taken to ensure that the cables are not kinked during sterilization as this can produce cracks in the cable and reduce cable life.

**ARCHIVED COPY:
DO NOT USE**

Cleaning the Generator

The TORS Generator may be cleaned as follows:

- i. Dilute a neutral pH detergent according to the manufacturer's directions.
- ii. Using above solution, lightly moisten a soft, clean cloth. Wipe surfaces of the generator Generator.
- iii. Using tap water, lightly moisten a soft, clean cloth. Wipe surfaces of the generator Generator.
- iv. Dry generator Generator surfaces with a soft clean cloth.

Cleaning the Footswitch

The TORS Footswitches may be cleaned as follows:

- i. Dilute a neutral pH detergent according to the manufacturer's directions.
- ii. Using above solution, lightly moisten a soft, clean cloth. Wipe surfaces of the footswitch Footswitch.
- iii. Using tap water, lightly moisten a soft, clean cloth. Wipe surfaces of the footswitch Footswitch.
- iv. Dry footswitch Footswitch surfaces with a soft clean cloth.
- v. Do not allow any water to enter into the air hoses.
- vi. Do not detach the air hoses from the footswitch Footswitch.

**ARCHIVED COPY.
DO NOT USE**

Formatted: Normal

Equipment Care

Probe Inspection

Before use it is recommended that the Probe be inspected for damage.

TORS Probes are susceptible to damage if forced into contact with metal (e.g. hand instruments and clamps, etc.) when active. The consequence of scratching a Probe is to raise the mechanical stress in the region of the scratch (creating a "stress raiser") when the Probe is vibrating. If the stress is raised sufficiently there is a risk of the Probe suffering metal fatigue and cracking. Fatigue failure is more likely if the stress raiser is close to one of the fixed points of maximum mechanical stress in the vibrating Probe. If a Probe does suffer metal fatigue as a result of a stress raiser causing a crack, the Probe will no longer possess a resonant frequency "recognizable" by the generator. Report to the supplier any Transducer with a Probe that has a scratch, if the scratch is deep enough to be detected by sliding a fingernail over it.

Formatted: Font: Italic

Formatted: Font: Italic

Formatted: Justified

Transducer / Probe Usage

The Cement and Soft tissue Transducers plus the titanium Probes of TORS can be decontaminated and reused. This is an important factor in reducing the cost per case. However, they cannot be reused indefinitely and a lifetime is set within which optimum performance can be expected. In order to monitor the so-called "lifetime" of the Transducer and Probe a memory chip is housed inside the cable plug of every Transducer.

For the Soft tissue transducer, every time the generator energizes the Transducer, the duration of the operation is monitored by the TORS then written to the memory chip, overwriting the previous total. When the total run time reaches pre-set markers for warning or termination, the TORS will alert the user.

For the Cement transducer, each time the Transducer is used for a procedure, the remaining number of allowed uses is monitored by the TORS then written to the memory chip, overwriting the previous total. When the total run reaches zero TORS will alert the user and further use will be prevented.

ARCHIVED COPY!

Storage of Equipment between Cases

It is recommended that the Transducers (with associated cables) are stored, between cases, in a large autoclave tray provided and again, care should be taken to ensure that the cables are not kinked close to the connector.

DO NOT USE

Maintenance

The hospital is responsible for ensuring that the unit has an electrical safety check performed by qualified service personnel at least once a year.

Do not remove the covers from TORS.

TORS generator does not require periodic calibration. If the Generator detects an internal problem it will display a "Service Due" on the rear LCD. If this is seen, contact RSL to arrange repair.

There are no user-serviceable parts in TORS.

Any damage to the Transducers or cables should be reported and the components returned to the supplier at the earliest opportunity.

Technical Support

In line with the Company's policy of quality assurance and customer care, RSL Ltd. offers a telephone helpline service for users of TORS. This is provided by the TORS manufacturer, RSL Ltd at its facility in Devon, England.

HELPLINE Telephone No. as displayed on the unit is: +44 (0)-1364- 653899

APPENDICES

Appendix 1: Markings on the TORS front and back panels

	<u>Cement transducer selected</u>
	<u>Soft tissue transducer selected</u>
	<u>Cement transducer output - Channel 1</u>
	<u>Cement transducer output - Channel 2</u>
	<u>Soft tissue transducer output - Channel 3 (low power)</u>
	<u>Soft tissue transducer output - Channel 3 (high power)</u>
	<u>Transducer active</u>
	Model <u>n</u> Number
	Serial <u>n</u> Number
	Manufacturer
	Date of <u>m</u> Manufacture
	<u>FUSE rating - Time delay, 5A, High breaking capacity, 250 Volts AC, size 20mm</u>
100-240V	Mains AC <u>v</u> oltage <u>range</u>
	Alternating <u>c</u> urrent
50/60Hz	Mains AC <u>f</u> requency
200VA	Input <u>p</u> ower
	Warning sign

ARCHIVED COPY:
DO NOT USE

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt









Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt

Formatted: Font: (Default) Arial, 9 pt




Formatted Table

	Type BF e Equipment
36kHz	Cement f Frequency of o Output
<70W	Cement o Output p Power
For Intermittent Use OFF/ON/OFF 20/20s	Cement d Duty c Cycle
36kHz	Soft t Tissue f Frequency of o Output
<50W	Soft t Tissue o Output p Power
For Intermittent Use OFF/ON/OFF 3/30s	Soft t Tissue d Duty c Cycle
	Volume
	Footswitch connection
	Equipotential c Connection
	Output channel
	Follow instructions for use
	Electrical and Electronic equipment. Return waste to a collection system or treatment and recycling facilities. Follow decontamination instructions before returning waste.
	General caution sign

ARCHIVED COPY:

DO NOT USE

Formatted Table

	Follow instructions for use
	Electrical and Electronic equipment. Return waste to a collection system or treatment and recycling facilities. Follow decontamination instructions before returning waste.
	General Caution Sign

ARCHIVED COPY:
DO NOT USE

Appendix 2: TORS (Series 1) Parts List

ITEM	Product Code
Cement Axial Grip Handpiece (-Single-Use-)	T1AH
Cement 200 Ø9 Fan Probe (-Reusable-)	F9R2T1FM-9-200
Cement 200 Ø8 Two-Lobe Groover Probe (Reusable) Cement 300 Ø9 Fan Probe (-Reusable-)	G8R2T1FM-9-300
Cement 200 Ø9 Extractor Probe (Reusable) Cement 200 Ø15 Fan Probe (-Reusable-)	E9R2T1FM-15-200
Cement 200 Ø6 Piercer Probe (Reusable) Cement 300 Ø15 Fan Probe (-Reusable-)	P6R2T1FM-15-300
Cement 200 Ø8 Piercer Probe (Reusable) Cement 200 Ø20 Fan Probe (-Reusable-)	P8R2T1FM-20-200
Cement 200 Ø10 Piercer Probe (Reusable) Cement 300 Ø20 Fan Probe (-Reusable-)	P10R2T1FM-20-300
Cement 200 Ø12 Piercer Probe (Reusable) Cement 200 Ø8 Scraper Probe (-Reusable-)	P12R2T1SM-8.5-200
Cement 300 Ø8 Scraper Probe (Reusable)	T1SM-8.5-300
Cement 300 Ø6 Piercer Probe (Reusable)	T1PM-6.5-200
Cement 300 Ø6 Piercer Probe (Reusable)	T1PM-6.5-300
Cement 200 Ø8 Piercer Probe (Reusable)	T1PM-8.5-200
Cement 200 Ø8 Scraper Probe (Reusable) Cement 300 Ø8 Piercer Probe (-Reusable-)	S8R2T1PM-8.5-300
Cement 200 Ø10 Scraper Probe (Reusable) Cement 200 Ø8 Piercer/Scraper Probe (-Reusable-)	S10R2T1P/SM-8.5-200
Cement 200 Ø12 Scraper Probe (Reusable) Cement 300 Ø8 Piercer/Scraper Probe (-Reusable-)	S12R2T1P/SM-8.5-300
Cement 200 Ø8 Two-Lobe Groover Probe (Reusable)	T1GM-8.5-200
Cement 300 Ø8 Two-Lobe Groover Probe (Reusable)	T1GM-8.5-300
Cement 200 Ø9 Extractor Probe (Reusable)	T1EM-9.5-200
Cement 300 Ø9 Extractor Probe (Reusable)	T1EM-9.5-300
Soft Tissue Transducer + Probe (-Reusable-)	T1TT
Cement Transducer (-Reusable-)	T1CT
Cement Cable (-Reusable-)	T1CC
Generator	T1G
Generator Cart	T1C
Mains Cable	234
Generator Carry Case	T1GC
Cement Footswitch	T1FC
Soft Tissue Footswitch	T1FT
Mains Cable	234
Autoclave Tray	T1AT
Autoclave Tray Mat	T1SM
Instrument Carry Case	T1IC

ARCHIVED COPY.
DO NOT USE

Spanner	T1S
Probe Cleaning Tool	T1CL
Probe Cleaning Cell	T1CW

[Optional Extras](#)

Formatted: Space Before: 18 pt, After: 6 pt

Cement 200 Ø14 Piercer Probe (Reusable)	P14R2
Cement 200 Ø16 Piercer Probe (Reusable)	P16R2
Cement 200 Ø18 Piercer Probe (Reusable)	P18R2
Cement 200 Ø20 Piercer Probe (Reusable)	P20R2
Cement 200 Ø16 Scraper Probe (Reusable)	S16R2
Generator Cart	T1C

ARCHIVED COPY:
DO NOT USE

Appendix 3: Electromagnetic Interference

Table 1

Guidance and manufacturer's declaration – electromagnetic emissions			
TORS T1GLG4 is intended for use in the electromagnetic environment specified below. The customer or the user of TORS should assure that it is used in such an environment.			
Emissions Test	Compliance	Electromagnetic environment - guidance	
RF emissions CISPR 11	Group 1	TORS uses RF energy only for its internal function. Therefore its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. TORS is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.	
RF emissions CISPR 11	Class A		
Harmonic emissions IEC 61000-3-2	No testing – not connection to public mains network		
Voltage fluctuations /flicker emissions IEC 61000-3-3	No testing – not connection to public mains network		

Table 2

Guidance and manufacturer's declaration – electromagnetic immunity			
TORS is intended for use in the electromagnetic environment specified below. The customer or the user of TORS should assure that it is used in such an environment.			
IMMUNITY test	EN 60721 test level	Comparison level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±6 kV contact ±8 kV air (1)	Floors should be conductive. No synthetic material should be used in the environment. The relative humidity should be in the range 40% to 60%.
Electrical fast transient/burst IEC61000-4-4	±2 kV for power supply lines	±2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	±1 kV differential mode ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % U_T (2) (100 % dip in U_T) for 0.5 cycle at: 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315° 0 % U_T (100 % dip in U_T) for 1 cycles 70 % U_T (30 % dip in U_T) For 25/30 cycles 0 % U_T (100 % interrupt in U_T) for 250/300 cycles	0 % U_T (100 % dip in U_T) for 0.5 cycle at: 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315° 0 % U_T (100 % dip in U_T) for 1 cycles 70 % U_T (30 % dip in U_T) For 25/30 cycles 0 % U_T (100 % interrupt in U_T) for 250/300 cycles	Mains power quality should be that of a typical commercial or hospital environment. If the user of TORS requires continued operation during power mains interruptions, it is recommended that TORS be powered from an uninterruptible power supply.
Power frequency (50/60 Hz) magnetic field IEC61000-4-8	3 A/m	No testing	No magnetically sensitive components.
<p>NOTES:</p> <p>Mitigation applied because of environment.</p> <p>U_T is the a.c. mains voltage prior to application of the test level.</p>			

Table 3


Guidance and manufacturer's declaration – electromagnetic immunity			
TORS is intended for use in the electromagnetic environment specified below. The customer or the user of TORS should assure that it is used in such an environment.			
IMMUNITY test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3V rms 150kHz to 80 MHz Outside ISM bands 60 rms ISM bands 0.15MHz to 80MHz 80% AM at 1kHz	3V rms	Portable and mobile RF communications equipment should be used no closer to any part of TORS, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.2\sqrt{P}$ $d = 1.2\sqrt{P}$, 80MHz to 800MHz $d = 2.3\sqrt{P}$, 800MHz to 2.3GHz where P is the maximum power output rating of the transmitter in Watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.3GHz 9 V/m 710MHz, 745MHz, 780MHz, 5240MHz, 5500MHz, 5785MHz	3 V/m 9V/m	Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. ^a
Immunity to proximity fields from RF wireless communications equipment	27 V/m 385MHz 28 V/m 450MHz, 810MHz, 870MHz, 930MHz, 1720MHz, 1845MHz, 1970MHz, 2450MHz	27V/m 28V/m	Interference may occur in the vicinity of equipment marked with the following symbol. 
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
<i>a Field strengths from fixed transmitters-, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the in the location in which TORS is used exceed the applicable RF compliance level above, TORS should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating TORS.</i>			
<i>b Over the frequency range 150 kHz to 80 MHz-, field strengths should be less than 3 V/m</i>			

Table 4

Recommended separation distances between portable and mobile RF communications equipment and TORS			
TORS is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of TORS can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment and TORS as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter	Separation distance according to frequency of transmitter		
	m		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800MHz to 2.5 GHz $d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance (in metres (m)) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

ARCHIVED COPY

DO NOT USE

Appendix 4: Tones and Beeps

Rising flourish	TORS is booting to standby mode
No Tone	TORS is in standby mode
Triple tone – low-medium-low pitch	Handpiece a Activate b Button has been pressed to initialize handpiece (Cement) Handpiece f Toggle b Button has been pressed to initialize handpiece (Soft Tissue)
Continuous low pitched tone	Acoustic output at the LOW power level
Continuous high pitched tone	Acoustic output at the HIGH power level
Triple tone – low-medium-high pitch	Soft f issue f Transducer is changing from LOW to HIGH power
Triple tone – high-medium-low pitch	Soft f issue f Transducer is changing from HIGH to LOW power
Double beep –high>low pitch	1. The G enerator has reset itself after a minor problem such as a mistuning event, time-out or over-temperature warning or 2. The G enerator has encountered a minor problem that requires power switch to be cycled OFF/ON
Triple beep – high>high>high pitch	A more serious problem has been detected

**ARCHIVED COPY:
DO NOT USE**

Appendix 5: Rear Panel Messages

The LCD on the rear panel will display messages that indicate the status of the equipment. The following table shows the possible displays:

Message	Status	Action required
Radley Scientific Ltd TORS Issue x	Start up message. Shows the Issue number of the software "X"	Now awaiting Transducer to be plugged in
Cement / Soft Tissue Ready	The Toggle switch has been operated and the relevant Transducer is ready to be activated.	Now awaiting use
36500Hz Low Power	Low power is selected (STT Only) (via Soft tissue footswitch or handpiece)	No action required
36500Hz High Power	High power is selected (STT Only) (via Soft tissue footswitch or handpiece)	No action required
36000Hz Cement Active OR Soft Tissue Active	While Transducer is active the top line will show the frequency. After releasing the switch it will display the final running frequency.	No action required
Active Too long Release Switches	The active button has been held on for too long. No output.	Release activate button on handpiece or on footswitchFootswitch
Release Switches	Either a Activate or f Toggle button (or F ootswitch es) has been pressed at switch on.	Release a Activate or f Toggle button when switching on Generator
Check Transducer Ease Grip and Heavy	Transducer has been loaded too heavily Transducer is too hot.	Release switch, then reactivate using less pressure on the jaw Allow transducer to cool.
Transducer Limit	Transducer is approaching 100% lifetime. Final use.	Replace Transducer after use
Change Tdcr Transducer Limit Transducer Limit Ready	Transducer has reached 100% lifetime and must be replaced.	Replace Transducer
Change Transducer Restart	Transducer frequency too low and feedback signal is low.	Switch supply off & on. Replace Transducer if seen 3 times
Transducer Leakage Change Tdcr	Generator has detected voltage on the Transducer.	Switch supply off & on. Replace Transducer if seen 3 times
Frequency Error Service due	The G enerator has detected a serious internal problem.	Switch supply off & on. If message seen again Generator requires service.

|

**ARCHIVED COPY:
DO NOT USE**

Appendix 6: Technical Specification

Model Nos:	See Appendix 2	
Dimensions:	Generator:	340mm (width) x 95mm (height) x375mm (depth) 13.4" (width) x 3.7" (height) x 13.4" (depth)
Weight:	Generator:	7.6 kg
	Transport case:	13.8 kg (loaded with generatorGenerator)
	Transducer:	0.37 kg
Fuse Type:	T5A, 250V, 20mm (2 off)	
Cordset	Contact RSL for recommended type	
Power supply input	100V - 240V, 50/60Hz	
Power consumption	200VA	
Cement Output - Frequency of operation	36kHz	
Cement Output – Power	<70W	
Cement Mode of Operation	Intermittent ON/OFF, 20/20s	
Cement Transducer Classification:	Type: BF	
Soft Tissue Output - Frequency of operation	36kHz	
Soft Tissue Output – Power	<50W	
Soft Tissue Mode of Operation	Intermittent ON/OFF, 3/20s	
Soft Tissue Transducer Classification:	Type: BF	
Accuracy of frequency display	1%	
Insulation Classification	Generator:	Class 1
Transducer	Titanium, s stainless steel, and plastic Autoclavable maximum 50 cycles/	
Environment for Transportation & Storage:	Temperature: -10°C to +50°C Relative humidity: 10% to 90% Atmospheric pressure: 50kPa to 106kPa	
Environment for Use:	Temperature: +10°C to +30°C Relative humidity: 30% to 75% Atmospheric pressure: 81kPa to 106kPa	

RSL Ltd. will make available on request circuit diagrams, component part lists, descriptions, calibration instructions to assist to service personnel in parts repair.

TORS has been designed and built in accordance with ISO 13485: 2016 Quality Assurance standard for medical devices.

CE conformance has been certified and the equipment complies with:

IEC 60601-1:2005 + CORR. 1:2006 + CORR.2:2007 + A1:2012
EN 60601-1:2006 + A11:2011 + A1:2013
ANSI/AAMI ES60601-1:2005/ (R) 2012

Appendix 7: Warranty Statement

Subject to the terms and conditions listed below, Radley Scientific Ltd. (hereafter called "the Company") guarantees to replace or repair free of charge any defective parts of TORS notified within the warranty period. This applies to the hardware defined below for the purposes of warranty claims made by any party supplied directly by Company.

Conditions of Warranty:

1. The Warranty covers the electronic system and ancillary components including the ~~generator~~Generator and ~~footswitch~~Footswitch, and applies to defects or faults that are reported to the Company, or its approved service facility, within ONE YEAR from the date of purchase.
2. Where transducers, cables or reusable probes are purchased separately or as part of the electronic system, the Warranty applies within the criteria of Clauses 5 and 6.
3. Where transducers, cables or reusable probes remain the property of the manufacturer the Warranty does not apply. If the Company loans a handset for use with the rest of the system the Company retains responsibility for performance.
4. The Warranty specifically excludes non-disposable (Handpieces and ~~Single Use~~ probes). The disposables are degradable consumables for use with the ~~Transducer~~ Transducer. They are for ~~Single Use~~ Single Use and are disposed of thereafter. Its condition and sterility is guaranteed up to the expiry date shown on the packaging label.
5. The Warranty relates to any components in Clauses 1 and 2 which are faulty or below specification as a result of inferior workmanship or materials. At the Company's expense the Company will conduct analysis, repair such faults, update to the relevant specification, re-test, re-calibrate and certify as necessary. For handsets this will take into account the additional criteria of Clause 6.
6. Transducer warranty criteria are defined as follows.
 - Cable, connector, acoustic system within (including waveguide) and transducer case are guaranteed against defective parts, materials or workmanship for 9 months from invoice date whether used or not, or until the transducer reaches end of life ("transducer limit") through use.
7. This Warranty does not cover accidental damage to, or fair wear and tear on, any part of the system.
8. This Warranty is only valid for use in the purchasing hospital and when used by personnel trained by the Company or its authorised distributors.
9. If a repair of a fault is required following expiry of the warranty, the Company will conduct analysis and submit a quotation to the customer stating recommended action and cost to the customer.
10. In the unlikely instance of the equipment requiring recall the Company will contact the purchasing hospital immediately and will, at the Company's expense, collect the relevant components.

|

**ARCHIVED COPY:
DO NOT USE**



ARCHIVED COPY:
DO NOT USE



Bremridge House,

Bremridge,

Ashburton

S. Devon TQ13 7JX

UK

Tel: (44) (1364) 653899 – Helpline

www.radleyscientific.co.uk

Page 55 of 56

**ARCHIVED COPY:
DO NOT USE**

